

**BROOKFIELD DIGITAL RHEOMETER**

**MODEL DV-III+    Operating Instructions**

**Manual No. M/98-211-A0701**

Please record the Model and Serial Number of your viscometer.  
Having this information readily available will help us to assist you  
should there be any questions regarding your instrument.

Model No. \_\_\_\_\_

Serial No. \_\_\_\_\_



SPECIALISTS IN THE  
MEASUREMENT AND  
CONTROL OF VISCOSITY

**BROOKFIELD ENGINEERING LABORATORIES, INC.**  
11 Commerce Boulevard, Middleboro, MA 02346-1031 USA

TEL 508-946-6200 or **800-628-8139** (USA excluding MA)  
FAX 508-946-6262 INTERNET [www.brookfieldengineering.com](http://www.brookfieldengineering.com)



# ***LAB Online Exhibition***



# Contents

<b>I. INTRODUCTION .....</b>	<b>3</b>
I.1 Components .....	4
I.2 Utilities .....	5
I.3 Specifications .....	5
I.4 Data Retention .....	5
I.5 Set-Up .....	6
I.6 Connections .....	9
I.7 Key Functions .....	10
<b>II. GETTING STARTED .....</b>	<b>12</b>
II.1 Autozero .....	12
II.2 Rheometer Display .....	13
II.3 Spindle Entry .....	15
II.4 Direct Speed Entry .....	16
II.5 External Control .....	17
<b>III. MAKING VISCOSITY MEASUREMENTS .....</b>	<b>19</b>
III.1 Quick Start .....	19
III.2 Preparation .....	19
III.3 Selecting a Spindle/Speed .....	20
III.4 Multiple Data Points .....	20
III.5 Cleaning .....	20
<b>IV. PROGRAMMING THE DV-III+ AND ANALYSIS .....</b>	<b>21</b>
IV.1 Programming Concept .....	21
IV.2 DV-III Speed/Time Pair Programming .....	22
IV.3 Bevis Programs .....	34
IV.4 Choosing the Best Data Collection Method .....	36
IV.5 Data Analysis .....	37
<b>V. OPTIONS .....</b>	<b>40</b>
V.1 Set Up .....	40
V.2 Print .....	41
V.3 Alarms .....	42
V.4 Set Temperature .....	43
V.5 Data .....	43
V.6 Timed data collection .....	44
<b>VI. RHEOLOADER USER'S MANUAL .....</b>	<b>47</b>
VI.1 Description of B.E.V.I.S. Commands .....	47
VI.2 Example Programs .....	50
<b>APPENDIX A - Cone/Plate Rheometer Set-Up .....</b>	<b>52</b>
<b>APPENDIX B - Viscosity Ranges .....</b>	<b>56</b>
<b>APPENDIX C - Variables in Viscosity Measurements .....</b>	<b>61</b>
<b>APPENDIX D - Spindle and Model Codes .....</b>	<b>63</b>
<b>APPENDIX E - Calibration Procedures .....</b>	<b>67</b>
<b>APPENDIX F - VS-27Y Clamp Assembly .....</b>	<b>74</b>
<b>APPENDIX G - DV-III+ to Computer Command Set .....</b>	<b>75</b>
<b>APPENDIX H - Fault Diagnosis and Troubleshooting .....</b>	<b>77</b>
<b>APPENDIX I - Warranty Repair and Service .....</b>	<b>80</b>

## I. INTRODUCTION

The Brookfield DV-III+ Programmable Rheometer measures fluid parameters of Shear Stress and Viscosity at given Shear Rates. Viscosity is a measure of a fluid's resistance to flow. You will find a detailed description of the mathematics of viscosity in the Brookfield publication "*More Solutions to Sticky Problems*", a copy of which was included with your DV-III+.

The principle of operation of the DV-III+ is to drive a spindle (which is immersed in the test fluid) through a calibrated spring. The viscous drag of the fluid against the spindle is measured by the spring deflection. Spring deflection is measured with a rotary transducer. The measuring range of a DV-III+ (in centipoise) is determined by the rotational speed of the spindle, the size and shape of the spindle, the container the spindle is rotating in, and the full scale torque of the calibrated spring.

There are four basic spring torque series offered by Brookfield:

<u>Model</u>	<u>dyne•cm</u>	<u>mN•m</u>
LVDV-III+	673.7	0.0673
RVDV-III+	7,187.0	0.7187
HADV-III+	14,374.0	1.4374
HBDV-III+	57,496.0	5.7496

The higher the torque calibration, the higher the measurement range. The measurement range for each torque calibration may be found in Appendix - B.

**All units of measurement are displayed according to either the CGS system or the SI system.**

1. Viscosity appears in units of centipoise (shown as "cP") or milliPascal-seconds (shown as mPa•s).
2. Shear Stress appears in units of dynes/square centimeter ("D/cm<sup>2</sup>") or Newtons/square meter ("N/m<sup>2</sup>").
3. Shear Rate appears in units of reciprocal seconds ("1/SEC").
4. Torque appears in units of dyne-centimeters or Newton-meters (shown as percent "%" in both cases).

The equivalent units of measurement in the SI system are calculated using the following conversions:

	<u>SI</u>	=	<u>CGS</u>
Viscosity:	1 mPa•s	=	1 cP
Shear Stress:	1 Newton/m <sup>2</sup>	=	10 dyne/cm <sup>2</sup>
Torque:	1 N•m	=	10 <sup>7</sup> dyne•cm

References to viscosity throughout this manual are done in CGS units.

**WARNING:** Use of this instrument in a manner not specified by Brookfield may result in incorrect readings or instrument failure. Please read this manual prior to using the instrument.

## I.1 Components

Component	Part Number
DV-III+ Rheometer	depends on model
Powerbase	DVP-2Y
includes:	
Leveling Screws (3)	VS-3
Upright Rod	VS-20
Jam Nut	VS-21
Clamp Assembly	VS-27Y
Spindle Set with Case	
<b>LVDV-III+</b> set of four spindles or	SSL
<b>RVDV-III+</b> set of seven spindles or	SSR
<b>HA/HBDV-III+</b> set of seven spindles	SSH
<i>For Cone/Plate versions: a spindle wrench, one cone spindle and sample cup Part No. CPE-44Y replace the spindle set.</i>	
Power Cord	
for 115 VAC	DVP-65
for 230 VAC	DVP-66
RTD Temperature Probe	DVP-94Y
Ribbon Cable	DVP-145
Guard Leg:	
<b>LVDV-III+</b>	B-20Y
<b>RVDV-III+</b>	B-21Y
Carrying Case	DVP-71Y
RHEOLOADER Software	DVP-201Y
Cable (DV-III+ to Computer)	DVP-80
Operator Manual	M/98-211
<i>Please check to be sure that you have received all components, and that there is no damage. If you are missing any parts, please notify Brookfield Engineering or your local Brookfield agent immediately. Any shipping damage must be reported to the carrier.</i>	

## I.2 Utilities

Autosensing Power Supply:

Input Voltage: 90 - 260 VAC  
Input Frequency: 50 - 60 Hz  
Power Consumption: Less than 220 UA

Power Cord Color Code:

	United States	Outside United States
Hot (live)	Black	Brown
Neutral	White	Blue
Ground (earth)	Green	Green/Yellow

## I.3 Specifications

Speed Range: 0.01-250 RPM, 0.01 RPM increments from 0.01 to 0.99 RPM, 0.1 RPM increments from 1.0 to 250 RPM

Viscosity Accuracy:  $\pm 1.0\%$  of full scale range for a specific spindle running at a specific speed.

Temperature sensing range: - 100°C to 300°C (-148°F to 572°F)

Temperature accuracy:  $\pm 1.0^\circ\text{C}$  from -100°C to 150°C  
 $\pm 2.0^\circ\text{C}$  from +150°C to 300°C

Analog Torque Output: 0 - 1 Volt DC (0 - 100% torque)

Analog Temperature Output: 0 - 4 Volts DC (10mv / °C)

Printer Output: Centronics, serial

Computer Interface: RS232

Weight: Gross Weight: 35 lbs. 15.9 kg  
Net Weight: 32 lbs. 14.5 kg


Carton Volume: 2.0 cu. ft. 0.057 m<sup>3</sup>

## I.4 Data Retention

The DV-III+ will save spindle parameters (used to calculate centipoise, shear rate and shear stress), default settings and the test data from the last program test run when the rheometer is turned off or there is a power failure.


## I.5 Set-Up

- 1) Place the upright rod into the hole at the front of the base. The rack gear and clamp assembly should face the rear of the base (see Figure 1). The upright rod is held in place with the jam nut which is attached from the bottom of the base. Tighten this nut with a suitable wrench (spanner).

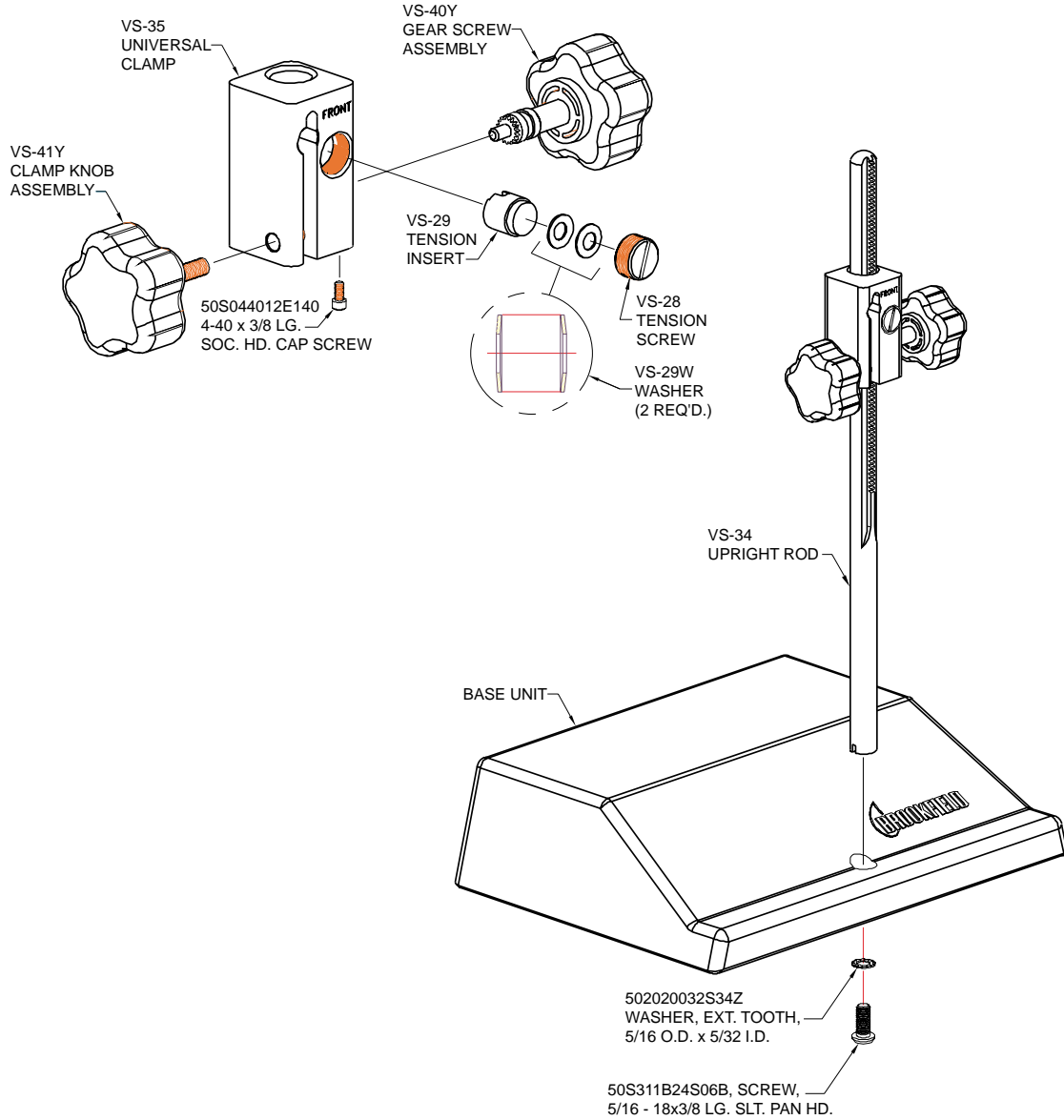

**Figure 1**

- 2) Insert the mounting handle on the back of the DV-III+ into the hole on the clamp assembly (Figure 2).

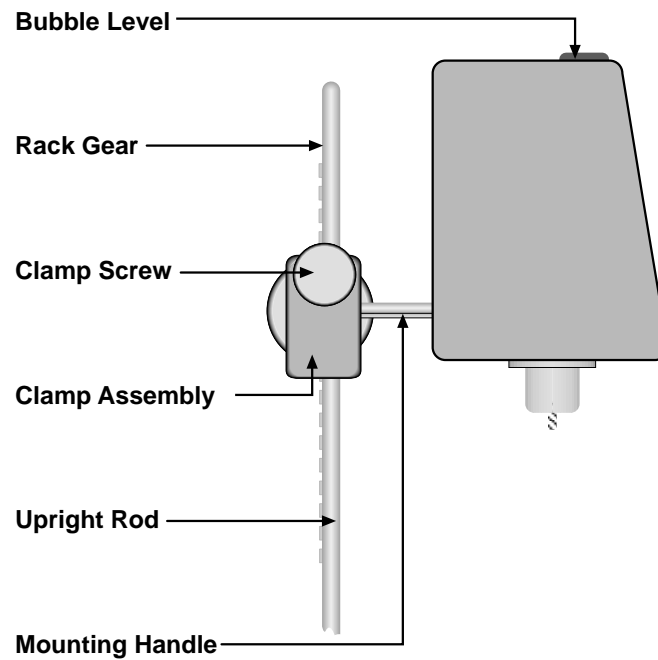


**Figure 2**

- 3) Tighten the DV-III+ clamp Screw (Figure 2).

**Note:** If the clamp assembly moves along the upright rod too freely, tighten the tension screw (see Appendix F).

- 4) Insert the ribbon cable into the DV-III+ Rheometer head. Insert the other end of the ribbon cable into the connector on the DV-III+ base (see Figure 3).

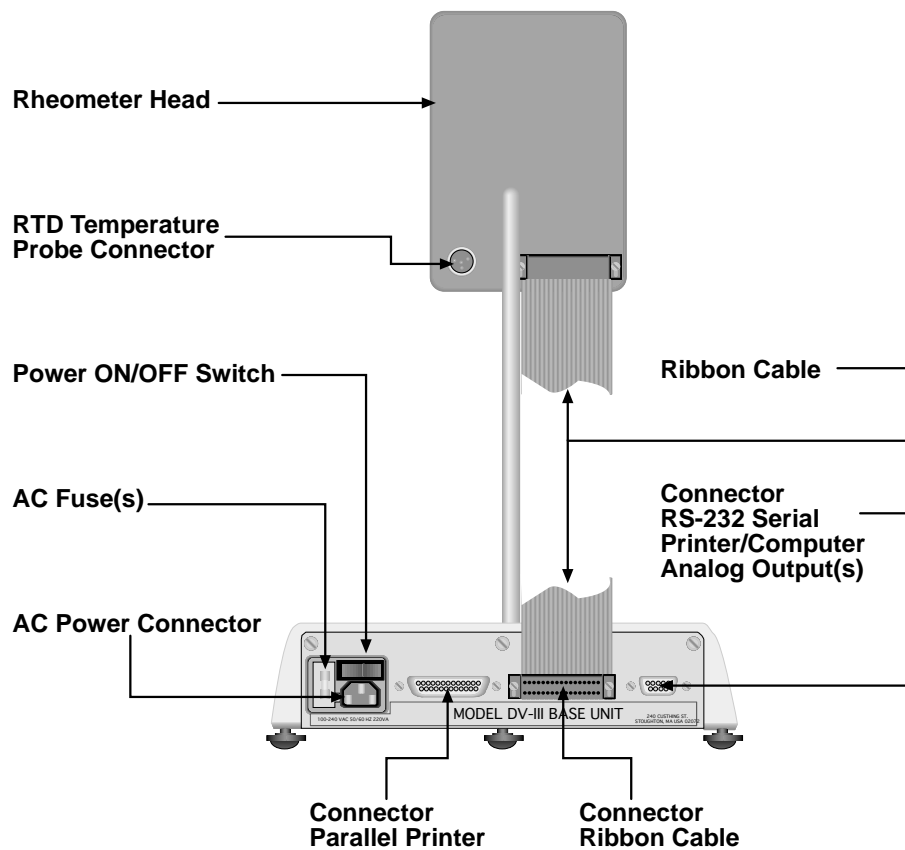

**Figure 3**

- 5) Connect the RTD probe to the socket on the back side of the DV-III+ Rheometer (Figure 3).
- 6) The Rheometer must be leveled before the instrument is zeroed and readings are taken. The level is adjusted using the three leveling screws on the base. Adjust so that the bubble level on top of the DV-III+ (Figure 2) is centered within the circle.
- 7) Make sure that the AC power switch at the rear of the base unit is in the OFF position. Connect the AC plug to the socket on the back of the DV-III+ base and plug it into the appropriate AC line.

**The DV-III+ must be earth grounded to ensure against electronic failure!!**

- 8) Temperature monitoring is assured (after the instrument has stabilized) to within  $\pm 1.0^{\circ}\text{C}$  in the range  $-100^{\circ}\text{C}$  to  $+150^{\circ}\text{C}$  and within  $2^{\circ}\text{C}$  in the range  $150^{\circ}\text{C}$  to  $300^{\circ}\text{C}$ .
- 9) For Cone/Plate models refer to Appendix A.
- 10) For printers, software and temperature controllers, refer to Section 1.6, Connections.

## I.6 Connections

The DV-III+ Rheometer is capable of communicating with several external devices to enhance operation. The cables and connections required for proper communication are detailed below.

### ◆ RHEOLOADER SOFTWARE

DVP-80 cable is used to connect the RS232 serial port on the DV-III+ base to Com Port 1 or Com Port 2 on the computer. This cable is supplied with the DV-III+.

### ◆ RHEOCALC SOFTWARE

DVP-80 cable is used to connect the RS232 serial port on the DV-III+ base to Com Port 1 or Com Port 2 on the computer. This cable is supplied with the RHEOCALC software.

### ◆ PARALLEL PRINTER

CAP-86 cable is used to connect the 25-pin parallel port on the DV-III+ base with the Centronics port on the printer.

### ◆ SERIAL PRINTER

DVP-81 cable is used to connect the 9-pin serial port on the DV-III+ with the 25-pin serial port on a printer.

### ◆ THERMOSEL CONTROLLER, MODEL HT-106

TC-200/TC-500/TC-201P/TC-501P BATH, MODEL HT-107

DVP-141 cable is used to connect the serial port on the DV-III+ base to the serial port on the controller. This cable is supplied with the controller/bath.

Be sure that the controller temperature probe is properly located in the control device (Thermosel or bath) and connected to the controller.

- Notes:**
1. The controller may alternately communicate with Rheocalc V 2.0 software. In this configuration, the controller is connected to the computer through either Com Port 1 or Com Port 2. The DV-III+ is also connected to a computer Com Port.
  2. The controller must also be connected to the control device (Thermosel or bath) with the appropriate load cable.

### ◆ STRIP CHART RECORDER

DVP-96Y cable is used to connect the serial port on the DV-III+ to the input block of the strip chart recorder. This cable is supplied with a Brookfield strip chart recorder.

## I.7 Key Functions

Figure 4 shows the control keys on the face of the DV-III+ Rheometer. The following describes each key's function.

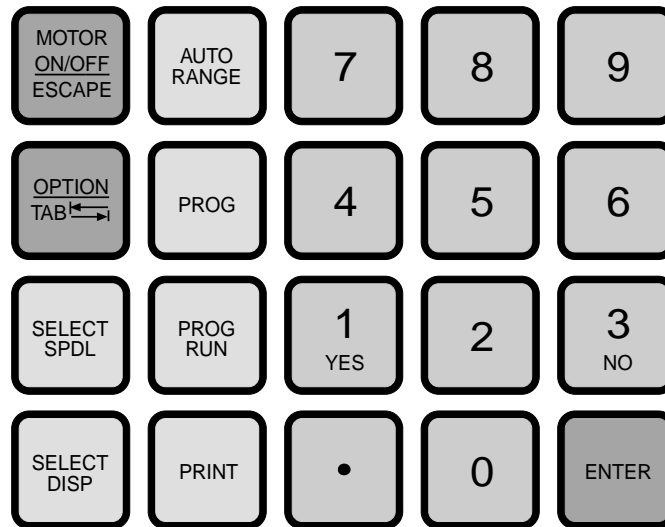


Figure 4



### MOTOR ON/OFF, ESCAPE

Turns the motor on or off. Cancels any operation, returns the user to the previous screen.



### AUTORANGE

Presents the maximum (100% torque) viscosity attainable using the selected spindle at the current speed.



### SELECT SPDL

Allows selection of the spindle to be used.



### SELECT DISP

Selects the parameter to be displayed:

- %** Rheometer Torque (%)
- cP** Viscosity (cP or mPa·s)
- SS** Shear Stress (Dynes/cm<sup>2</sup> or Newtons/m<sup>2</sup>)
- SR** Shear Rate (1/Sec)



### OPTION, TAB

Accesses options menu. See Section V. Toggles between selectable items when indicated.



### PRINT

Sends a single line of data to an attached printer. Selects printing and non-printing mode as selected in the Options menu.



## PROG

Access the Programs menu for program creation, running or deleting. Constructs a test program. Allows you to review/modify an existing test program. Execute a Bevis program.



## PROG RUN

Execute DV-III speed/time pair program.



## NUMBER KEYS (0 through 9)

Sets speeds and choose items from various dialog screens and the option menu.



## ENTER

Functions as an ENTER key similar to a computer by serving to accept a keyboard entry.

## II. GETTING STARTED

### II.1 Autozero

Before readings may be taken, the Rheometer must be autozeroed. This is done each time the power switch is turned on. The Rheometer will guide you through the procedure, as follows:

Turn power switch on; as shown in Figure 5, the screen indicates that the DV-III+ is in the standalone mode (is not connected to a computer) and gives the version of the operating firmware (the built in program which controls the instrument) and a two-digit alphanumeric code which indicates the Model number (see Table D2 in Appendix D; the code tells the spring torque rating of your Rheometer).



Figure 5

No key press is necessary. After a short pause the display will read “*REMOVE SPINDLE, LEVEL RHEOMETER AND PRESS THE MOTOR ON/OFF KEY TO AUTOZERO.*” Before beginning the autozero procedure, Brookfield recommends that you allow 10 minutes for the instrument to warm up.

After pressing the MOTOR ON/OFF key, the screen “flashes” for approximately 15 seconds while the DV-III+ autozeros.

After 15 seconds the display reads “*AUTOZERO IS COMPLETE REPLACE SPINDLE AND PRESS ANY KEY.*” Press a key.

The main screen is displayed and the DV-III+ is ready for use (Figure 6).

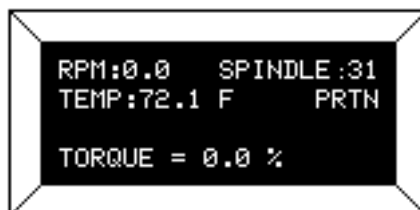


Figure 6

## II.2 Rheometer Display

The DV-III+ Rheometer is supplied with a 4-line display. The basic set of information is called "The Default Screen" and is shown in FIGURE 7. The parameters are detailed below:

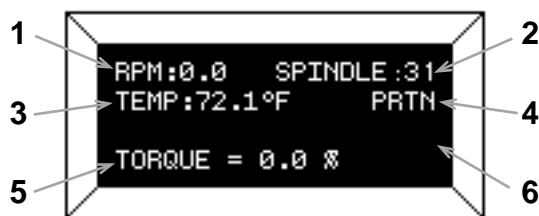


Figure 7

### 1. Motor Status and Current Rheometer Speed

The DV-III+ motor can be OFF, ON at 0.0 rpm or ON at a speed greater than 0.0 rpm. When the motor is OFF, "OFF" will be displayed and no speed entry will be accepted. When the motor is ON, the actual speed of rotation will be displayed. When the motor is switched from ON to OFF, the speed of rotation will be remembered; when the motor is turned ON again, the DV-III+ will operate at that same speed. The rheometer motor is set to "OFF" after AUTOZERO.

**Note:** Motor OFF and a speed setting of 0.0 are essentially the same.

### 2. Spindle Number

The currently-selected spindle. Viscosity, shear rate, and shear stress values will be calculated based on this number. See Section II.3.

### 3. Measured Temperature

The current temperature as measured by the attached temperature probe. If no probe is connected, four dashes "----" will be displayed.

### 4. Printing Status

Indicates the currently-selected method of printing. See Section II.5.

### 5. Measured Data

Instrument Torque (%), Viscosity (cP), Shear Stress (D/cm<sup>2</sup>), Shear Rate (s<sup>-1</sup>)  
The parameters are toggled from one to another using the Select Display key.

**Note:** Shear Stress and Shear Rate data cannot be calculated for some spindle geometries. In these cases, the display will show 0.0.

### 6. Blank Line

This line is used to display entry data when selecting a spindle or speed of rotation. Additionally, selected programs available for running will be identified here when in the Program mode. (See Section IV.2).

The default screen will appear at the completion of the AUTOZERO sequence each time the DV-III+ is turned ON in the standalone mode (see Section II.6 for external control mode). The displayed data may be changed as described in the following sections.

The format for data displayed in the default screen and all other screens is described in Table 1. For appearance sake, the entries in the table have been decimal point aligned. Actual rheometer display will have all fields left justified.

Item	Print	Format	Range	Example
RPM	RPM	X.XX XX.X XXX.X	0.01 <= RPM <= 0.99 0.1 <= RPM <= 99.9 100 <= RPM <= 250	0.09 2.4 150.0
Model	M	XXXXXX	See Model Table D2 - App D	RV
Spindle	S	XX	00 <= S <= 99	31
Torque	T	XX.X	-10.1 <= T <= 99.9	82.4
Viscosity	cP or mPas	X.XX XXX.X XXXXXX XXXeX	0 <= cP <= 9.99 10 <= cP <= 999.9 1000 <= cP <= 99999 100000 <= cP <= 51200000000	3.16 123.8 12345 123e3 to 8
Shear Stress	D/CM2 or N/M2 = D / C M 2 / 1 0	X.XX XXX.X XXXXXX XXXeX	0 <= D/CM2 <= 9.99 10 <= D/CM2 <= 999.9 1000 <= D/CM2 <= 99999 100000 <= D/CM2 <= 999999	4.56 234.5 12345 123e3
Shear Rate	1/SEC	X.XXX XX.X XXXXXX	0 <= 1/SEC <= 9.999 10 <= 1/SEC <= 99.9 100 <= 1/SEC <= 99999	1.234 20.7 200
Temperature	T		-99.9 <= T <= 300.0	-10.3
Time	Z	XX:XX	00:00 <= Z <= 99.59	05:32

**Table 1**


◆ OUT OF RANGE INDICATORS

The DV-III+ is capable of measuring instrument torque within the range of 0 to 100%. Based on this measurement, viscosity and shear stress are calculated. Brookfield recommends that data be collected only in the range of 10 to 100%. Any data collected outside of this range is considered invalid.

The DV-III+ provides the following display indicators when the measurement point is outside of the 10-100% acceptable range.

◆ TORQUE GREATER THAN 100%

When Rheometer torque exceeds 100%, the parameter display field will show “EEEE” for torque, viscosity and shear stress.

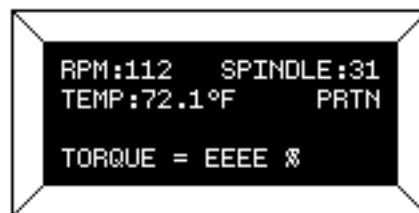

**Figure 8**

◆ TORQUE LESS THAN 10%

When Rheometer torque drops below ten (10) percent, the Rheometer will *continue* to display measurement (% , cP, D/cm<sup>2</sup>) values with units flashing:

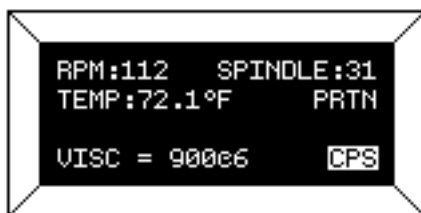


Figure 9

◆ TORQUE LESS THAN 0%

When Rheometer torque drops below zero (0) percent, the Rheometer will *continue* to display torque values preceded by a minus (-) sign. The viscosity and shear stress field will display dashes (- - - -) as indicated in the next screen display:

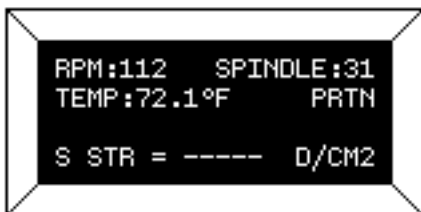


Figure 10

### II.3 Spindle Entry

The user can elect to change the spindle selection by pressing the SELECT SPDL key. The DV-III+ control program will use the previously blank line 3 on the default display screen to record the new spindle input as depicted in **Figure 11**.

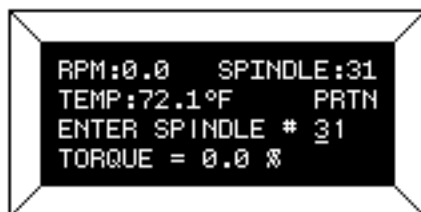


Figure 11


To enter a spindle number, press the numeric keys until the desired spindle number has been entered. Valid spindle numbers encompass the range from 00 to 99 as listed in **Appendix D**. Mistakes are corrected by repeatedly pressing the numeric keys until the proper spindle value has been entered. At that point, the user presses the SELECT SPDL key again. An invalid spindle entry will result in a “beep” and the display of the data entry error screen as depicted below.



Figure 12

An invalid spindle entry is any two digit number in the range from 01 to 99 which is not listed in **Appendix D**. This error message will be displayed for a few seconds after which the spindle entry screen (**Figure 11**) will be re-displayed with a blank field for the spindle number. The user can cancel spindle entry at any time by pressing the MOTOR ON/OFF/ESCAPE key.

The user may elect to use a special spindle whose selection is accomplished by first entering a spindle number of **99** and then pressing the SELECT SPDL key. This will result in the following display:

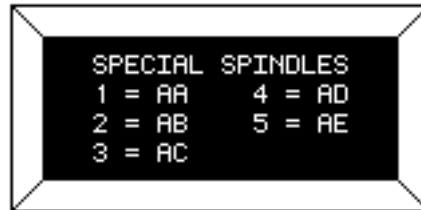

**Figure 13**

At this point, press the numeric key for the special spindle of choice. This list is created at the time the Rheometer is manufactured. This list will therefore depend on the number of special spindles ordered and could contain as few as one (1) or as many as five (5) spindles. If no special spindles were purchased, the following message will be displayed if **99** is entered for a spindle number:




**Figure 14**

Press any key to exit this screen and to return to the spindle selection screen. The user may again select another spindle or press the SELECT SPDL key to cancel spindle selection operation.

Successful selection of a spindle at the press of the SELECT SPDL key returns the user to the default screen with the new spindle displayed in the upper right-hand corner. For standard spindles this would be the two (2) digit designator used to select the spindle. In the case of special spindles, the two (2) letters (**AA, AB, AC, AD** or **AE**) corresponding to the special spindle would be displayed instead. The spindle number or letters will be retained in memory when power is removed. This means that the last value entered for the spindle will be displayed the next time the Rheometer is turned on.

## **II.4 Direct Speed Entry**


At this point, the user may choose to enter a speed by the so-called direct speed entry method. Enter a valid speed in the range of **0.01** to **250 RPM** by pressing the numeric keys successively. The previously blank line 3 on the default display screen records the user's new speed input as depicted in Figure 15:

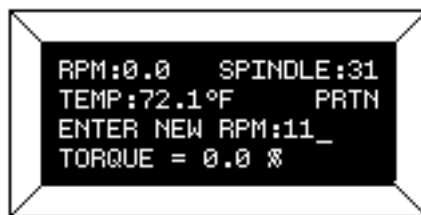

**Figure 15**

Here, the user intends to enter a speed of **112 RPM**, has pressed the “1” key twice and is about to press the “2” key. If the user makes more than five (5) key presses, the **DV-III+** control program will “roll” the cursor back to the first character of the field and begin to overwrite the previous data entry.

Next the user presses the ENTER key to accept the speed. The motor will begin running at **112 RPM** and the display will be updated to the next screen image:

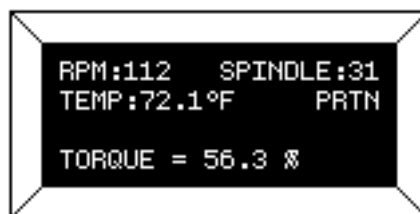

**Figure 16**

If the speed entered was not valid the Rheometer will display the following message:



**Figure 17**

After a few seconds, the display returns to **Figure 15** with the speed data field cleared and just the underscore cursor awaiting a new entry.

## II.5 External Control

The DV-III+ Rheometer can be used in conjunction with Brookfield software, RHEOCALC (V. 2. or higher). Through RHEOCALC, all rheometer functions are controlled by the computer. The DV-III+ must be set to the external control mode to allow for proper communication with RHEOCALC. To configure the external control mode, connect cable DVP-80 to the serial port on the DV-III+ base **before** turning on the DV-III+. With the DVP-80 cable in place, the DV-III+ will present the screen shown in Figure 18 when it is turned on. If external control is selected, the DV-III+ will display Figure 19 and only accept control commands from RHEOCALC software.



**Figure 18**



**Figure 19**

The DV-III+ may be set to stand alone mode by turning it OFF and ON again and selecting "Stand Alone" or by removing the DVP-80 cable prior to turning the DV-III+ on.

**Note:** The DV-III+ cannot communicate with RHEOLOADER software in the external control mode. Chose "Stand Alone" when presented with Figure 18 if you want to use RHEOLOADER.

## III. MAKING VISCOSITY MEASUREMENTS

### III.1 Quick Start

The DV-III+ Rheometer uses the same methodology as the Brookfield Dial Reading Viscometer and DV series of Digital Viscometers. If you have experience with other Brookfield equipment, this section will give you the quick steps for taking a viscosity reading. If you have not used a Brookfield Viscometer before, skip this section and go to Section III.2 for a detailed description.

- A) Assemble and level the rheometer (Section I.5).
- B) Autozero the rheometer (Section II.1).
- C) Enter the spindle number using the SELECT SPINDLE key (Section II.3).
- D) Introduce the spindle into the sample and attach the spindle to the coupling nut.  
**NOTE: Left-hand threads.**
- E) Enter the speed of rotation using the number pad and ENTER key (Section 11.4).
- F) Record % torque and viscosity.

### III.2 Preparation

- A) **RHEOMETER:** The DV-III+ should be turned on, leveled and autozeroed. The level is adjusted using the three feet on the bottom of the base and confirmed using the bubble on the top of the head. Adjust the feet until bubble is inside the center target. Set the level prior to autozero and check the level prior to each measurement.

**Proper level is essential for correct operation of the DV-III+.**

- B) **SAMPLE:** The fluid to be measured (sample) must be in some container. Many spindle systems from Brookfield are supplied with specific sample chambers such as the Small Sample Adapter, UL Adapter and Thermosel. The standard spindles supplied with the DV-III+, LV (1-4), RV (1-7) and HA/HB (1-7), are designed to be used with a 600ml low form Griffin beaker (or equivalent container with a diameter of 8.25 cm).

Brookfield recommends that you use the appropriate container for the selected spindle. You may choose to use an alternate container for convenience, however, this may have an effect on the measured viscosity. The DV-III+ is calibrated considering the specified container. Alternate containers will provide results that are repeatable but not "true."

The LV (1-4) and RV (1-7) are designed to be used with the guardleg attached. Measurements made without the guardleg will provide repeatable results but may not provide "true" results.

**When comparing data with others, be sure to specify the sample container and presence/absence of the guardleg.**

Many samples must be controlled to a specific temperature for viscosity measurement. When conditioning a sample for temperature, be sure to temperature control the container and spindle as well as the sample.

Please see our publication, "More Solutions to Sticky Problems", for more detail relating to sample preparation.

### III.3 Selecting a Spindle/Speed

The DV-III+ has the capability of measuring viscosity over an extremely wide range (for example, the RVDV-III+ can measure fluids within the range of 100-40,000,000 cP) (see Appendix B). This range is achieved through the use of several spindles over many speeds.

The process of selecting a spindle and speed for an unknown fluid is normally trial and error. **An appropriate selection will result in measurements made between 10-100 on the instrument % torque scale.** Two general rules will help in the trial and error process.

- 1) Viscosity range is inversely proportional to the size of the spindle.
- 2) Viscosity range is inversely proportional to the rotational speed.

In other words: to measure high viscosity, choose a small spindle and/or a slow speed. If the chosen spindle/speed results in a reading above 100%, then reduce the speed or choose a smaller spindle.

Experimentation may reveal that several spindle/speed combinations will produce satisfactory results between 10-100%. When this circumstance occurs, any of the spindles may be selected.

Non-Newtonian fluid behavior can result in the measured viscosity changing if the spindle and/or speed is changed. See our publication, "More Solutions to Sticky Problems," for more detail.

**When viscosity data must be compared, be sure to use the same spindle, speed, container and temperature.**

### III.4 Multiple Data Points

The majority of viscosity measurements are made at the quality control level and consist of a single data point. The test is conducted with one spindle at one speed. The data point is a useful bench mark for the go/no-go decision in a production setting. The DV-III+ can be used for single point measurement.

Many fluids exhibit a characteristic change in viscosity with a change in applied force. This non-Newtonian flow behavior is commonly seen in paints, coatings and food products as a decrease in viscosity as shear rate increases. This behavior cannot be detected or evaluated with the single viscosity point measurement.

Non-Newtonian flow is analyzed through the collection of viscosity data over a range of shear rates and the generation of a graph of viscosity versus shear rate (a rheogram). This information will allow for a more complete characterization of a fluid and may help in formulating and production of the product. The DV-III+ is capable of collecting multiple data points for the analysis of flow behavior. See Section IV on Programming and Analysis.

More information on flow behavior, shear rate and rheograms is available in our publication, "More Solutions to Sticky Problems."

### III.5 Cleaning

All immersed components are stainless steel. Use cleaning solutions that are not corrosive and avoid scratching the measurement surfaces. The instrument housing should be cleaned with a soft damp cloth.

## IV. PROGRAMMING THE DV-III+ AND ANALYSIS

The programming and data analysis functions of the DV-III+ are accessed by pressing the PROG key on the rheometer. The display will change to present a menu with three choices: DV-III, B.E.V.I.S., and Models. DV-III and B.E.V.I.S. are the programming alternatives. Models will present the five math models available for data analysis.

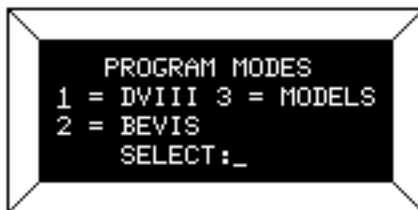


Figure 20

### IV.1 Programming Concept

The DV-III+ may be programmed to collect viscosity data without operator involvement. The captured data may be displayed and analyzed or output to a printer. Programs may be written using two different methodologies, DV-III and B.E.V.I.S.

The DV-III programming technique uses speed/time pairs to control the DV-III+. A program consists of multiple lines (up to 25) instructing the rheometer to operate at a particular speed for some period of time. As an example, we can instruct the DV-III+ to rotate the spindle at 5 RPM for 30 seconds and then change speed to 10 RPM and wait 20 seconds with the following program:

Step 1	RPM = 5	Time = 00:30
Step 2	RPM = 10	Time = 00:20

A single data point will be collected at the end of each time interval.

The B.E.V.I.S. programming technique uses a custom program language to control the DV-III+. A program consists of a series of commands instructing the rheometer in speed control, time control, data collection, temperature control, and output. B.E.V.I.S. offers a higher level of rheometer control compared to the DV-III method. However, the construction of B.E.V.I.S. programs is more involved. The 2-step DV-III program previously described is duplicated using B.E.V.I.S. commands below:

```
SSN    5
WTI    00:30
PDN
SSN    10
WTI    00:20
PDN
END
```

The involved programming of B.E.V.I.S. commands is a small trade for the significant increase in control capability over the DV-III method.

## IV.2 DV-III Speed/Time Pair Programming

This programming method allows the operator to control the DV-III+ through the variables of speed and time. These speed/time pairs instruct the rheometer to operate at a speed of rotation for a certain period of time. Programs can be created with up to 25 steps. The DV-III+ can store up to 10 programs. Upon completion of a program, the data may be viewed on the DV-III+ display, analyzed or printed to an attached parallel or serial printer.

Two examples of programs are shown below:

Collect Data Over Time			Collect Data At Several Speeds		
Step	RPM	Time	Step	RPM	Time
1	100	00:12	1	2.5	01:00
2	100	00:12	2	5	00:30
3	100	00:12	3	10	00:30
4	100	00:12	4	20	00:15
5	100	00:12	5	50	00:15

Five viscosity data points will be collected over one minute.

Five viscosity data points will be collected at five speeds over 150 seconds.

This program mode is accessed by pressing the program key and selecting number 1; 1 = DV-III. The creation, editing and execution of DV-III programs are described in the following sections.

There are two types of test programs:


- 1) **Next Speed Set** where the test speeds are programmed, and the operator must signal the DV-III+ to change speeds (and therefore take a reading) by pressing the ENTER key.
- 2) **Prog Speed Set** where the DV-III+ will perform the test automatically.

Each step of a program has two variables - speed and hold time. The reading is taken at the end of the hold time interval in a Prog Speed Set or when the ENTER key is pressed in a Next Speed Set.

If the first step hold time interval is 0 seconds, the program is a Next Speed Set type. If the first step interval is 1 second or more, the program is a Prog Speed Set type.

### ◆ SPEED SET SELECTION AND PROGRAMMING

The DV-III+ viscometer allows for the retention of a maximum of **10** speed sets with up to **25** discrete speeds per speed set. The program locations are numbered 0 through 9. These speed sets are retained in **EEPROM** memory for those times when the DV-III+ is not powered up. To access a previously programmed speed set or to enter data for a new speed set, the user presses the “1” key when in the display of **Figure 20** and is presented with the screen shown in **Figure 21**:

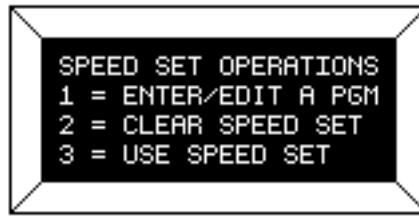

**Figure 21**

At this point, the user may Enter/Edit, Clear or Use a stored program (Speed Set). Let's start with Enter/Edit by pressing the "1" key:



**Figure 22**


In this example, the user is informed that he has 6 speed sets (0,1,2,5,8,9) pre-programmed in memory and 4 speed sets (3,4,6,7) not programmed and available. Select any one of the ten speed sets by pressing the appropriate numeric key. Pressing the **MOTOR ON/OFF/ESCAPE** key at this point would exit the user to the default **PROGRAM MODES** display (**Figure 20**). For now let's assume that the user wants to program a new speed set by pressing the "3" key (the first available program slot).

◆ ENTERING A SPEED SET (PROGRAM)

There are two (2) types of programs available to the user: programs with finite step time intervals and programs with zero (0) step time intervals. We will cover the inputting of finite step time programs first.

◆ SPEED SETS WITH FINITE STEP TIMES (PROG SPEED)

These programs when executed will automatically progress from step to step based on the time intervals programmed by the user. On pressing the "3" key in **Figure 22** the user is presented with:

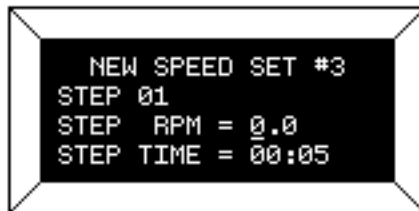

**Figure 23**

This screen reminds the user of the speed set that he has selected to program and then allows him to change either the speed or time interval or both for that step.

**Note:** The time interval on entry to this screen will always be set to 00:05 seconds as the default value. The user may of course change it to any valid time of his choice. Whenever you change time interval, that new time becomes the default interval until it is again changed by the user. Also, note that zero (0) times are not allowed for program steps after the first step for **Finite Step** programs.

The OPTION/TAB key is used for moving from input field to input field and the ENTER key to accept the current input for a step. On entry to this screen, the underscore cursor would be flashing (as shown) under the first digit of the step RPM. Use the numeric keys to make changes to the step speed, repeating the input as many times as required until satisfied.


When satisfied with the speed input, press the OPTION/TAB key which moves the flashing cursor down to the first character of the time field. The same procedure is used here to input the step time as was used to enter the RPM above. Speed or time data that is out of range, as defined by **Table 1**, will result in the following screen:



**Figure 24**


This screen will be displayed for **1-3** seconds.

When ready, the user may press the PROG RUN key to display the data for the next step in the program, or the MOTOR ON/OFF/ESCAPE key (whereby none of the changes up to that point will be accepted) to return to the screen of **Figure 22**. To end a program, the user simply enters and accepts a step RPM and Time of zero (0) or continues to input step data until the program reaches the twenty-five (25) step program limit. In either case, the following screen will be displayed:

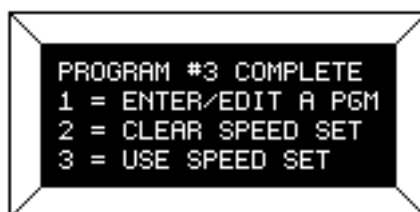

**Figure 25**

To use the currently selected speed set, press the “3” key in **Figure 25**. This would immediately revert to the default screen modified as follows:

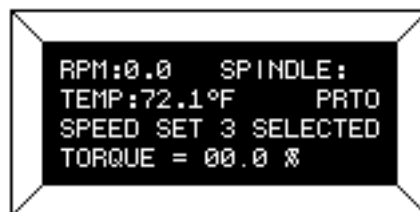


**Figure 26**

Note: If at this point, prior to using the program the user wished to enter a direct speed, a press of any numeric key which would result in a display similar to **Figure 15**. At the completion of the direct speed input, the display would revert to **Figure 26** above with the appropriate **RPM** displayed, and the viscometer running at that speed.

The program is initiated by pressing the PROG RUN key. See "Using Pre-programmed Speeds."

#### ◆ SPEED SETS WITH ZERO STEP TIMES

These programs when executed will require that the user press the ENTER key to progress from step to step. On pressing the "3" key in **Figure 22** the user is presented with the same screen that he saw in the above description for finite step programs:

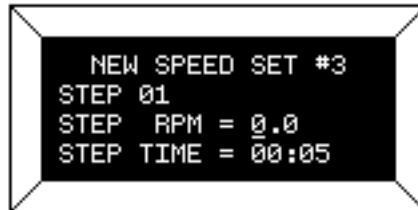


**Figure 27**

The user inputs his step RPM exactly as he did for finite step time programs above. However, for time, input 00:00 and press the ENTER key. From this point forward, the user will only be able to enter speeds since each press of the ENTER key will advance him to the next step. The OPTION/TAB key will not be required. If the user wishes to correct the speed input, continue to press the numeric/decimal point keys until satisfied. To correct a speed after pressing the ENTER key for that step, wait until the program is complete and then edit the program to correct the mistake. To end a program, simply enter and accept a step RPM of zero (0) or continue to input step data until the program reaches the twenty-five (25) step program limit. Speed restrictions/limits are the same as for the description just above as are the error messages.

#### ◆ EDITING A SPEED SET (PROGRAM)

This item is used to review a just-entered program or to review/modify (edit) a program already stored in a memory slot. Entry to this method would typically be from **Figure 22** after selecting an "IN MEM" program slot or by pressing the "1" key in **Figure 21** having just finished entering a program. In either case, the user is presented with:

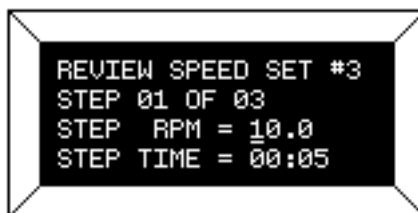

**Figure 28**

Operation in this mode is exactly the same as for entering a new speed set; all key actions and speed and time limits are the same. At this point, the user may continue to review/modify the speeds comprising speed set #3 or elect to print a listing of the speeds in this speed set. To accomplish this, the user must be in the program Enter/Edit mode; have selected or programmed a speed set which contains more than two (2) speeds, and then press the PRINT key. If all is well (i.e. satisfied the above requirements) the rheometer will display:



**Figure 29**

This message simply asks the user to make sure the printer is ready (it's on-line and has paper in it) and then awaits for the **PRINT** key to be pressed. When it is pressed, the **DV-III+** will send the following data to the attached printer:

◆ FOR SPEED SETS WITH FINITE STEP TIMES

Program Use: \_\_\_\_\_

Programmer: \_\_\_\_\_

BROOKFIELD DV-III+ RHEOMETER — DATA FOR SPEED SET #9

SPEED #01	RPM = 2.5	TIME INTERVAL = 00:05
SPEED #02	RPM = 5.0	TIME INTERVAL = 00:05
SPEED #03	RPM = 10.0	TIME INTERVAL = 00:05
SPEED #04	RPM = 20.0	TIME INTERVAL = 00:05
SPEED #05	RPM = 50.0	TIME INTERVAL = 00:05

◆ FOR SPEED SETS WITH ZERO STEP TIMES

Program Use: \_\_\_\_\_

Programmer: \_\_\_\_\_


BROOKFIELD DV-III+ RHEOMETER — DATA FOR SPEED SET #8

SPEED #01	RPM = 10.0	TIME INTERVAL = 00:00
SPEED #02	RPM = 20.0	TIME INTERVAL = 00:00
SPEED #03	RPM = 30.0	TIME INTERVAL = 00:00
SPEED #04	RPM = 40.0	TIME INTERVAL = 00:00

Since speed sets can contain twenty-five (25) separate speeds, printing the speeds that comprise a speed set will be of great help in allowing the user to fully exercise the power of the **DV-III+**. After the printing is complete, the user will be returned to the display of **Figure 28**.

◆ CLEARING A SPEED SET FROM MEMORY

Since **10** speed sets can be retained in memory, the user may eventually use all the available speed set slots. The user may also have programmed speed sets that are no longer required and would like to remove. Assume that the user had programmed a new speed set #3 above. If he wished to **permanently** remove that new speed set, or any other **IN MEMORY** speed set he would, while in the screens of **Figure 21**, press the “2” key and be presented with:

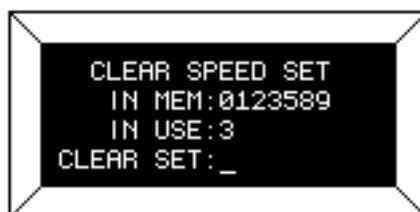

**Figure 30**


This screen advises that there are **7** speed sets in memory; speed set **#3** is in use and that the **DV-III+** is awaiting input for the speed set to delete.

Note: If no speed set is in use the word “**NONE**” will appear next to the **IN USE:** prompt.

At this point, the user has two options:


1. Pressing the **MOTOR ON/OFF/ESCAPE** key will exit from this screen and no speed sets will be cleared. *Or*
2. Pressing any of the keys “**0**”, “**1**”, “**2**”, “**3**”, “**5**”, “**8**”, or “**9**” will delete that speed set.

Thus, to discard speed set **#5**, the user would press the “**5**” key and be presented with:

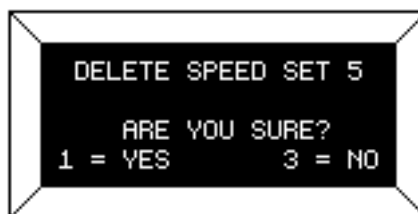

**Figure 31**

In which the **DV-III+** is requesting that the user specifically press the “**1**” key in order to delete the desired speed set.

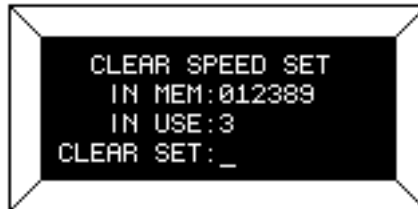

**Figure 32**

Pressing the “**3**” key will cause the **DV-III+** to take no action will return the user to the **CLEAR SPEED SET** opening screen, **Figure 30**. Any attempt to delete an in-use speed (“**3**” for instance) will cause the **DV-III+** to issue a “beep” “beep” with no action being taken. Thus no active (i.e. selected for use) program can be deleted from this screen.

#### ◆ USING PRE-PROGRAMMED SPEEDS

Pressing the three (**3**) key from **Figure 21** takes the user to the speed set selection screen of **Figure 22** where the user selects a new speed set. That done the user is sent to the default screen with: “**SPEED SET X SELECTED**” displayed on line three (3) of the screen.

The user initiates the use of programmed speeds by pressing the **PROG RUN** key. If the user presses the **PROG RUN** key with no speed set selected, the following error box will be displayed:



Figure 33

However, we will assume at this point that we have selected speed set #2 for use in the ensuing data gathering operations. To initiate the use of this speed set (with finite step times or with zero step times), the user presses the PROG RUN key and is presented with a start/end step input screen as shown next:

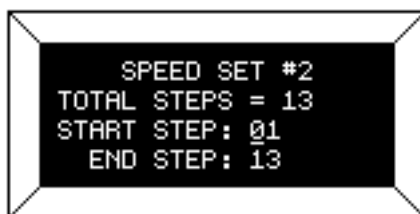


Figure 34

If the user had not previously entered start and end steps, this screen will display **01** for the start step, and the last program step (**13** in this case) as the end step. The user could elect to use the entire speed set at this point by pressing the PROG RUN key. If the user had previously selected a start and end step, those values would be displayed upon entry to this screen instead of the program limit values as shown above. However, while a speed set can contain up to **25** separate speeds, the user may be in a situation where only a few **contiguous** steps may be required. Therefore, this screen allows for the option of entering the range of speeds encompassed by the start step (not necessarily the first step) and the end step (not necessarily the last step). Pressing any numeric key at this point will erase the currently displayed start step and substitute the new value. The user may select a start step less than the end step ( a so-called **Up Ramp**) or a start step greater than the end step ( a so-called **Down Ramp**). Any attempt to enter a start or end speed not contained in the speed set will result in the following display:



Figure 35

The ENTER key is used to step from the start step entry to the end step entry. Repeated pressing of the ENTER key will allow the user to move back-and-forth between the start and end entries and change them as required until the correct start and end step values have been entered. The selected speed set, and the start and end steps values entered, will be retained in EEPROM memory for use the next time the viscometer is powered up in the stand-alone mode. Pressing the OPTION/TAB key at this point results in the following screen display:

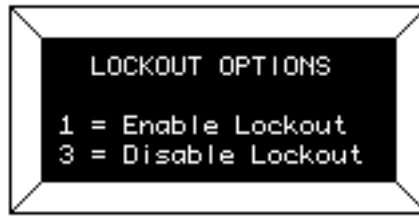


Figure 36

#### ◆ LOCKOUT OPTIONS

Pressing the “1” key **locks out** any use of the NUMERIC keys, and the PROGRAM, SELECT SPINDLE and OPTIONS/TAB keys. Pressing the “3” key would disable an existing lockout condition only when in the **LOCKOUT OPTIONS** screen, **Figure 36**. After pressing the “1” or “3” keys the user would be returned to the display of **Figure 34**. From **Figure 34** the user runs the program by pressing the PROGRAM RUN key which signifies that the user is satisfied with his start and stop step values, and wishes to start running with the selected speed set. Or, he can return to the default screen of **Figure 7** by pressing the MOTOR ON/OFF/ESCAPE key. (**Note:** the OPTIONS/TAB key is re-enabled at this point **only** to allow the user to return to the **LOCKOUT OPTIONS** screen). In either case, if the user has enabled the lockout mode, the top line of the default display will change as shown in **Figure 37** below:

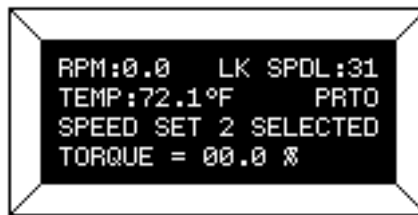


Figure 37

The “**LOCKOUT**” condition will remain in effect until the user re-starts the program mode by pressing the PROGRAM RUN key re-initiating the steps of **Figures 34** through **37** above.

Once the start and end steps have been selected, the program is started by pressing the **PROG RUN** key. However, if the viscometer motor was not turned on, the following message would be displayed:

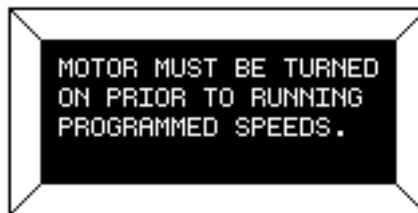


Figure 38


At this point, the user would turn the viscometer motor on by pressing the **MOTOR ON/OFF/ESCAPE** key. This will cause the viscometer to start running at the first selected program speed.

#### ◆ USING PROGRAMMED SPEEDS WITH A ZERO TIME INTERVAL

Assume the following:


- A speed set has been selected.
- A subset of the speed set has not been selected. (i.e. we will use the entire set of speeds)
- The speed set included a time interval that was equal to zero minutes and zero seconds.
- Print mode has been set to non-continuous mode (i.e. output will be sent to the printer only when the user presses the ENTER key).

With the above items in effect our default display screen will appear as follows after the **PROGRAM** steps have been completed:

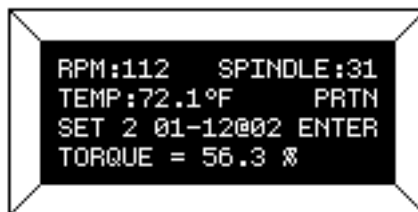

**Figure 39**

Or, if the user had selected the **LOCKOUT** option, the default screen above would appear as in **Figure 40** below:

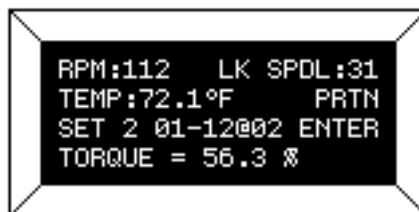


**Figure 40**

Here, the **LOCKOUT** mode is indicated by the revision to the top line of **Figure 40** where “**SPINDLE**” has been contracted to “**SPDL**” and the two-letter combo: “**LK**” is being used to indicate the “**LOCKOUT**” condition. Please note that while **TORQUE** is presently being displayed, there is nothing to prevent the user from pressing the SELECT DISP key to change the data display item.

The item to note here is the message being displayed on line 3 informing the user that the program is at speed #02 of the 12 speeds comprising speed set #02. The **DV-III+** will continue to operate at **112 RPM** until the ENTER key is pressed. At that point (pressing the ENTER key), two distinct events will occur:

1. The **DV-III+** will ramp up (or down) to programmed speed #03 and,
2. Since the printer is **On (PRTN)**, the data corresponding to the last speed (**112 RPM**) that existed at the time the ENTER key was pressed will be sent to the attached printer.

Thus, if 12 speeds comprised the selected speed set, then the user would be required to press the **ENTER** key 12 times to exercise the entire speed set. After the last speed has been executed (i.e. the user pressed the ENTER key for the twelfth time), the viscometer speed will be set at **ZERO RPM** and the following message will be displayed:

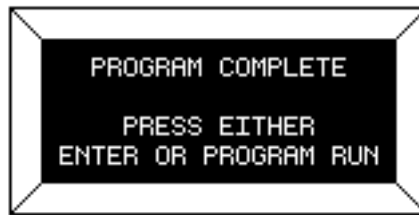

**Figure 41**


Pressing the ENTER key would return the user to the default screen, as depicted in **Figure 40**, or to the program start and stop limit selections of **Figure 34** if the PROG RUN key is pressed.

◆ USING PROGRAMMED SPEEDS WITH A NON-ZERO TIME INTERVAL

Assume the following:

- A speed set has been selected.
- A subset of the speed set has not been selected. (i.e. we'll use all of the speeds in the speed set)
- The speed set included a time interval that was greater than zero (0) minutes and zero (0) seconds.
- Print mode has been set to non-continuous mode (i.e. output will be sent to the printer only when the user presses the ENTER key).

With the above items in effect our default display screen will appear as follows after the program steps have been completed:

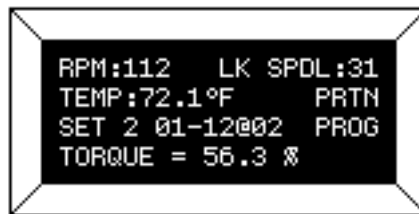


**Figure 42**

The difference between this mode and the zero time interval mode is reflected in the message being displayed on line three of the display where **PROG** has replaced **ENTER** indicating that the program will be executed step-by-step without user intervention. Here the user is again informed that he is at speed #02 of the 12 speeds comprising speed set #02. As the programmed time interval elapses, the following will occur:

- The **DV-III+** will ramp up (or down) to programmed speed #03 and,
- The data corresponding to the last speed (speed #02) that existed at the moment the step time elapsed will be sent to the printer since the printer is **ON (PRTN)** (The automatic ramping to the next speed will be interpreted by the **DV-III+** as if the ENTER key had been pressed, causing the printer output).
- At the completion of the speed set, the viscometer speed will be set at **ZERO RPM**

After the last speed has been executed, the viscometer speed will be set at **ZERO RPM** and the following message will be displayed:

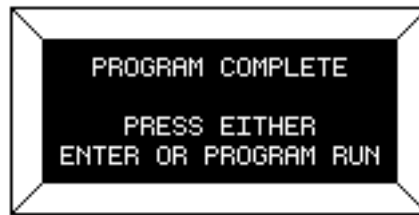


**Figure 43**

Pressing the ENTER key would return the user to the default screen, as depicted in its general form in **Figure 42**, or to the program start and stop limit selections of **Figure 34** if the PROG RUN key is pressed.


◆ OPERATION WITH PROGRAMMED SUBSET SPEEDS

Had the user selected a subset of a speed set (**Figure 34**), say speeds #04 through #08, then the programmed mode (with a non-zero time interval) would result in the following display:

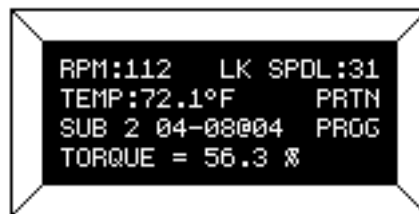


**Figure 44**

Here, the user is informed that speed set #2 is being used starting with speed number four and ending at speed #08 while currently executing speed #04. After the first speed (number #04) is completed, the display would be updated to show the new speed and the step display would now reflect the current executing step number as in **Figure 45** below.

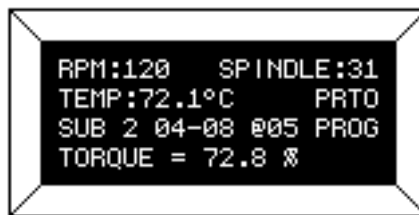

**Figure 45**

If the speed set contained **no** time interval, the selection of a subset of speeds would result in a screen display as shown in **Figure 46**.

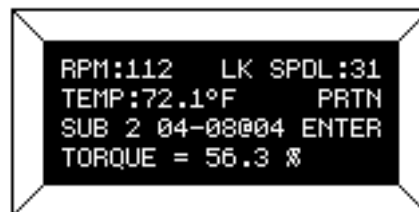

**Figure 46**

The step number would be incremented each time the user pressed the ENTER key. Assume that the speed corresponding to step #04 was executing. If the user presses the ENTER key, our display will be updated as shown in **Figure 47**



**Figure 47**

◆ PROGRAMMED SPEED STOP

The user may stop program mode operation at any time by pressing the MOTOR ON/OFF/ESCAPE key anytime during program operation. The following message would be displayed:




**Figure 48**

The viscometer motor would be automatically turned **OFF** and the default screen display (**Figure 37**) would show zero (0.0) **RPM**. At this point, the user may perform any valid viscometer operation - load a new speed set; run a direct speed; set alarms or even **re-start** the current program.

◆ PROGRAMMED SPEED HOLD

The user may interrupt program mode operation at any time by pressing the **PROG** key during program execution. The following message would be displayed:

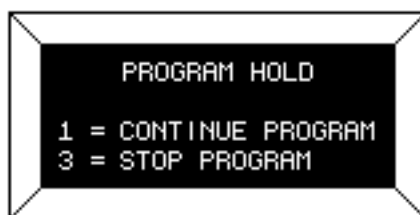

**Figure 49**

The viscometer motor would be automatically turned **OFF** and the **DV-III+** would be awaiting user key input. If the user presses the “1” key, the **DV-III+** will continue executing the current program, picking up from whichever **RPM** (and time, if applicable) it had been interrupted. However, if the user had pressed the “3” key, operation would be identical to that described above for **Programmed Speed Stop** mode.

### IV.3 Bevis Programs

The B.E.V.I.S. Programming Method allows the operator to control the DV-III+ through the variables of speed, temperature and time while providing for independent data collection. Programs can include up to 25 commands with a maximum data count of 800. The DV-III+ can store up to 10 programs. Upon completion of the program the data may be viewed on the DV-III+ display, analyzed or printed to an attached parallel or serial printer.

B.E.V.I.S. programs are created on a PC using Rheoloader software (supplied with the DV-III+). See Section VI for details. The programs are "loaded" onto the DV-III+. Loaded programs cannot be deleted, but can be overwritten.

The B.E.V.I.S. program menu is accessed by pressing the **PROG** key and selecting number 2, 2 = B.E.V.I.S. The loading and execution of B.E.V.I.S. programs are described in the following sections.

#### ◆ B.E.V.I.S. PROGRAMS MENU

**B.E.V.I.S.** operations are accessed by pressing the "2" key when in the **PROGRAM MODES** menu. The user is immediately presented with:

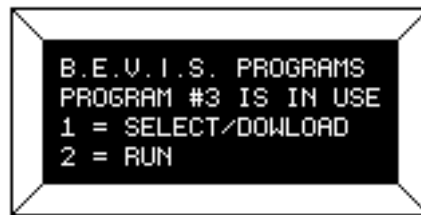


Figure 50

This screen informs the user that the **B.E.V.I.S.** program in storage slot **3** is current ("last used") and that it may be run by pressing the "2" key or another program may be selected by pressing the "2" key. It should be noted that entrance to the **B.E.V.I.S.** program mode makes the last used program available for printing or running. Thus the user, seeing that a program was resident in slot **3**, could have printed it directly from the above screen by pressing the front panel PRINT key.

#### ◆ SELECTING AND DOWNLOADING B.E.V.I.S. PROGRAMS

To select a program the user presses the "1" key and is presented with:

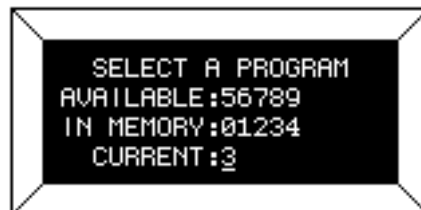


Figure 51

The user is informed that the current (or "last used") program is no. **3**; that programs **0,1,2,3** and **4** are available for immediate use and that five (5) slots: **5,6,7,8** and **9** are empty and are available for download from a host computer. A press of the ENTER key would select the current program (i.e. "3") for use while a press of an appropriate numeric key **0,1,2,3** and **4** would select that specific

program slot for use. The above screen would be updated to reflect the new selection. A press of the **5,6,7,8** or **9** keys will place the user in the **B.E.V.I.S.** program download mode. The following screen appears (assuming a press of the “**5**” key):



**Figure 52**


The user’s selection “**5**” is flashing and is the current slot selection. When ready, the user presses the ENTER key to begin the program download. The following screen will appear for the duration of the download. See Section VI. RHEOLOADER for information on creating and downloading B.E.V.I.S. programs.



**Figure 53**

◆ **RUNNING B.E.V.I.S. PROGRAMS**

Programs are run by pressing the “**2**” key when in **Figure 50** which presents the user with the following screen:

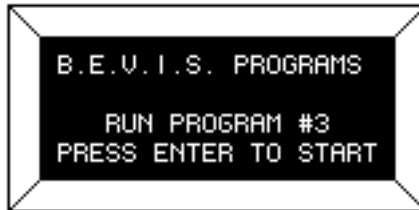

**Figure 54**


Here we see that our program is no. “**3**” and that a press of the ENTER key will start it running. As soon as the program starts executing the user will be presented with:



**Figure 55**

This is a typical display for a given program step. The current step being executed (**WTI0330**) is shown as well as the next step to be executed (**PDN**). The bottom line displays the current program

step, the total number of steps (**02/39**) and any time intervals if they are relevant. The **WKY** command message could be displayed on the bottom line in lieu of the step and time info. This screen stays resident until the user presses the **OPTIONS/TAB** key which “toggles” back-and-forth between this screen and an amended default screen shown next:

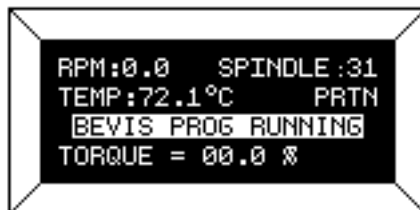

**Figure 56**

The user now sees viscosity data and can use the **SELECT DISPLAY** key to view other viscosity measurement parameters. Pressing the **OPTIONS/TAB** key from now on will toggle between the screens of **Figure 55** and **Figure 56**. The program code will automatically switch back to the program progress screen, **Figure 55** above, if a conditional has been reached, the end of a program step is reached or user input is required. If no user input is required, the **B.E.V.I.S.** program proceeds to the next step without switching back to the progress screen.

#### **IV.4 Choosing the Best Data Collection Method**

The DV-III+ offers 3 methods for data collection; Single Speed, DV-III Speed/Time Pairs and B.E.V.I.S. Programs. The decision of which technique is best should be made considering the test requirements.

##### ◆ SINGLE SPEED

Single speed measurements may be made by direct speed commands on the DV-III+ keypad. Viscosity and % Torque are read directly from the display. This technique offers the simplicity of the Brookfield Dial Viscometer. Multiple data points may be gathered by issuing multiple speed commands.

This technique is fast and easy. It's well suited to gathering data on samples prior to establishing a test method or for performing single point tests.

##### ◆ DV-III SPEED/TIME PAIRS

DV-III speed/time pairs offers a simple technique to collect multiple data points. Programs can be created and executed from the DV-III+ keypad. Results can be analyzed or output to a printer.

This technique is useful when multiple data points are required and the test method is simple.

##### ◆ B.E.V.I.S. PROGRAMS

B.E.V.I.S. Programs offer a command set capable of sophisticated rheometer DV-III and data collection. Programs are created on a PC and executed from the DV-III+ keypad. Results can be analyzed or output to a printer.

The B.E.V.I.S. Program technique is useful when sophisticated data collection is necessary.

## IV.5 Data Analysis

Data collected from DV-III speed/time pairs or B.E.V.I.S. programs may be analyzed using several math models. These models provide a means to numerically describe the behavior of the test fluid. In the case of viscosity measurement, a non-Newtonian fluid will produce a curve when test data is plotted on a shear stress vs. shear rate graph. The math model will force the data into a straight line and describe it with a slope and y intercept. The terminology associated with the slope and y intercept vary from model to model as does the interpretation of results.

The DV-III+ does not allow for data sets to be edited. Programs must be constructed to conform with the following data requirements if math models are to be used:

- The data set must contain non zero values for shear stress and shear rate (except for the paste model which requires non zero viscosity and RPM).
- There cannot be two equal adjacent shear rate values (RPM values for paste model).
- % torque values of all data points must be between 0.1% and 100%.

If any of the above circumstances are violated, an error message will appear when a math model is selected.


**ERROR #1:** A % torque value is less than 0.1. A shear stress or shear rate value is zero.

**ERROR #2:** A % torque value is greater than 100.

**ERROR #3:** Reserved


**ERROR #4:** Two adjacent speeds of equal value.

Math models for data analysis are accessed by pressing the PROG key and then 3; 3 = Models (Figure 20). If no data (i.e. no data at all or less than two(2) data points) is in the data buffer, no modeling can be performed and Figure 57 will be displayed:

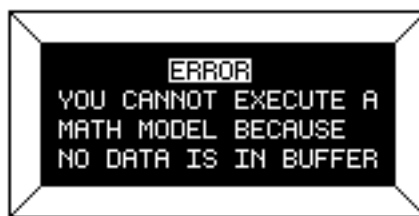


**Figure 57**

If there is already data in the data buffer, the user will be presented with the following screen:

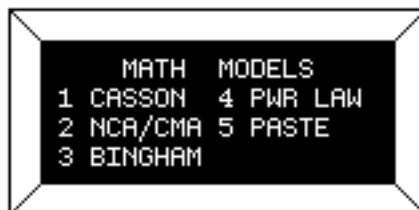


**Figure 58**

The user is informed that there are five (5) math models which can be used on the buffer data. A model is selected by pressing the appropriate numeric key. No matter the model selected, the following screen will be displayed for the duration of the mathematical analysis:

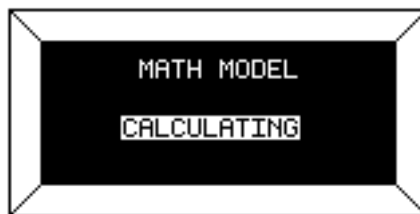

**Figure 59**

When the calculations are complete, the results for the particular model will be displayed as follows:

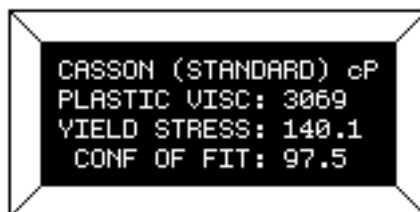

**Figure 60**

This screen, for the Standard CASSON Model, is typical for all five (5) of the math models. Note the cP in the upper right-hand corner to remind the user that the values are cP based. It could have been SI (if the user had opted for SI display in the **SETUP** menu) as shown on the next two (2) of the following four (4) math model screens:

The user may elect to print test results in order to obtain hard copy results. This is accomplished by pressing the **PRINT** key while any of the above screens are being displayed causing the following to be printed:

◆ STANDARD CASSON PRINTOUT

Sample Name: \_\_\_\_\_

Operator Name: \_\_\_\_\_

Date: 01/14/1999 Time: 02:27 Math Model Results: CASSON (STANDARD)

Model: HB Spindle: 34

Plastic Viscosity:	1906.3	cP
Yield Stress:	1976.88	D/CM2
Confidence of Fit:	63.6	%

The equation for each model is described below with a definition of parameters. Please contact Brookfield or an authorized representative if further information is required.

Note: The **confidence of fit** parameter used in all of the models is an indication of how well the model fits the data set. 100% indicates the best fit.

### 1. Casson (Standard)

The Standard Casson equation is:  $\sqrt{\tau} = \sqrt{\tau_0} + \sqrt{\eta D}$  where:

- $\tau_i$  = Shear Stress
- $\tau_0$  = Yield Stress (stress at zero shear rate)
- $\eta$  = Plastic Viscosity
- $D$  = Shear Rate

The calculated parameters for this model are:

- Plastic Viscosity (cP or mPa•s)
- Yield Stress (Dynes/cm<sup>2</sup> or N/m<sup>2</sup>)
- Confidence of Fit (%)

The Standard Casson method is a direct implementation of the original Casson equation.

### 2. NCA/CMA Casson (Chocolate)

This Casson method is derived from the standard set forth by the National Confectioners Association (NCA) and the Chocolate Manufacturers Association (CMA). Although based on the original Casson equation, this implementation has been tailored by the NCA and CMA specifically to applications involving chocolate.

The Chocolate Casson equation is:  $(1+a)\sqrt{\tau} = 2\sqrt{\tau_0} + (1+a)\sqrt{\eta D}$  where:

- $\tau_i$  = Shear Stress
- $\tau_0$  = Yield Stress (stress at zero shear rate)
- $\eta$  = Plastic Viscosity
- $D$  = Shear Rate
- $a$  = spindle (or bob) radius/ inner cup radius

The calculated parameters for this model are:

- Plastic Viscosity (cP or mPa•s)
- Yield Stress (Dynes/cm<sup>2</sup> or N/m<sup>2</sup>)
- Confidence of Fit (%)

### 3. Bingham Plastic

The Bingham equation is:  $\tau = \tau_0 + \eta D$  where:

- $\tau_i$  = Shear Stress
- $\tau_0$  = Yield Stress (stress at zero shear rate)
- $\eta$  = Plastic Viscosity
- $D$  = Shear Rate

The calculated parameters for this model are:

- Plastic Viscosity (cP or mPa•s)
- Yield Stress (Dynes/cm<sup>2</sup> or N/m<sup>2</sup>)
- Confidence of Fit (%)

### 4. Power Law

The Power Law equation is  $\tau = kD^n$  where:

- $\tau_i$  = Shear Stress
- $D$  = Yield Stress (stress at zero shear rate)
- $k$  = Plastic Viscosity
- $n$  = Shear Rate

The calculated parameters for this model are:  
Flow Index (no units)  
Consistency Index (cP or mPa•s)  
Confidence of Fit (%)

## 5. IPC Paste Analysis

This method is intended to calculate the Shear Sensitivity Factor and the 10 RPM Viscosity value of pastes. A prime example of its use is in the solder paste industry, thus the name IPC (Institute for Interconnecting and Packaging Electronic Circuits).

The Paste equation is:  $\eta = kR^n$  where:

$\eta$  = Viscosity (cP)  
 $k$  = Consistency Multiplier  
 $R$  = Rotational Speed (RPM)  
 $n$  = Shear Sensitivity Factor

The calculated parameters for this model are:

Shear Sensitivity Factor (no units)  
10 RPM Viscosity (cP or mPa•s)  
Confidence of Fit (%)

## V. OPTIONS

The DV-III+ Options Menu allows the user to execute temperature control commands and special time tests. General rheometer settings are also accessed from this menu. The Options Menu is shown in Figure 61. Selections are made by pressing the appropriate number key.



Figure 61

### V.1 Set Up

1. Change the units of temperature. The change is selected by pressing the 1 key. The change must be confirmed by pressing the ENTER key.
2. Change the units of viscosity and shear stress (CGS System: cP, D/cm<sup>2</sup>) (SI System: mPa•s, N/m<sup>2</sup>).

The change is selected by pressing the 2 key. The change must be confirmed by pressing the ENTER key.

3. Change communication status with external Brookfield temperature controller. An "off" indication means that there is no communication with a controller. Selecting this option will make the DV-III+ try to establish communication. When communication is established, "off" will be change to "on". Temperature control will always be set to "off" when the DV-III+ is turned on.

When communication is established, Line 1 of the default screen will be modified. The temperature field will show "CTLR" in place of "TEMP".

Communication may only be established with Brookfield controllers:

Thermosel Controllers:	HT-106, HT-104
Bath Controllers:	HT-107, HT-105

The change is selected by pressing the 3 key. Subsequent key presses required will be indicated on the display.

4. Change the temperature display by offsetting the measured temperature to agree with an external temperature measurement device. The adjustment will be indicated by flashing temperature units (F or C) on the default screen (**Figure 7**). The adjustment will be reset to 0.0 when the DV-III+ is turned on.

The adjustment is entered using the number keys. The sign (+ or -) is selected using the **OPTION/TAB** key. The change must be confirmed by pressing the **ENTER** key.

## V.2 Print

1. Change the time interval that is used when the DV-III+ is printing continuously. Data is entered in the format of MM:SS. For example: an interval of one minute and 30 seconds is entered as: 01:30.

The change is initiated by pressing the 1 key. Once the time interval is input, it must be accepted by pressing the ENTER key.

2. Change the port to be used for printing, parallel or serial.

The change is selected by first pressing the 2 key and then the appropriate key for the printing method. The change will be indicated by the position of the arrow on the right side of the display. The change is accepted by pressing the ENTER key.

When selecting serial printing, the DV-III+ will also ask for handshaking status. Please see the instruction manual of the printer for specification.

If the DV-III+ is communicating with an external temperature controller, serial printing will not be available(the controller uses the serial port).

3. Enter the date and time to be shown on printed data tables. Data is entered in the format of DD/MM/YY for date and HH:MM for time (24-hour clock). For example: 2:30 pm on January 15, 1999 is entered as 15/01/99, 14:30.

The change is initiated by pressing the 3 key. Once the data is input, it must be accepted by pressing the ENTER key.

## V.3 Alarms

There are three adjustable alarm settings: **LO ALARM %**, **HI ALARM %** and **MOTOR OFF %**. The values are set in the Set Alarms mode. Alarms are used to signal the operator that the fluid is out of the input specification. The alarms are set in % torque values, not Viscosity, Shear Stress or Shear Rate values. The range of values which may be entered for each alarm and their default values are:

LO ALARM	HI ALARM	MOTOR OFF
Minimum value: 10%	Minimum value: 0%	Minimum value: 0%
Maximum value: 99.9%	Maximum value: 100%	Maximum value: 115%
Default value: 10%	Default value: 99.9%	Default value: 110%

The procedure for entering and enabling alarm values is as follows starting from the main screen:

- 1) Press the OPTION/TAB key to display the Options Menu.
- 2) Press the NUMBER 3 key to display the Alarms Options screen.
- 3) Press the NUMBER 1 key to View/Set Alarms.
- 4) Enter the LO ALARM % torque value. The new entry will overwrite the default values.
- 5) Press the OPTION/TAB key to move the cursor to the HI ALARM % field.
- 6) Enter the HI ALARM % torque value. The new entry will overwrite the default values.
- 7) Press the OPTION/TAB key to move the cursor to the MOTOR OFF % field.
- 8) Enter the MOTOR OFF % torque value. The new entry will overwrite the default value.
- 9) Press the ENTER key to accept the ALARM values.
- 10) Press the NUMBER 3 key to enable (turn on) / disable (turn off) the alarms.
- 11) Press ENTER to accept the ALARM condition.

**Note:** The LO ALARM is tripped after the % torque reading falls below the setting. The beeping may be shut off by either the % torque reading rising above the alarm setting or by pressing the ESCAPE key.

The HI ALARM is tripped after the % torque reading goes above the alarm setting. The beeping may be shut off by either the % torque reading falling below the alarm setting or by pressing the ESCAPE key.

The MOTOR OFF is tripped after the % torque reading goes above the motor off setting. The DV-III+ stops rotating and the instrument beeps. Pressing any key turns off the beep.

## V.4 Set Temperature


The DV-III+ can issue temperature control commands when an external temperature control device has been connected (see Section V.1). Selecting item #4 Set Temp in the Options menu will display the current setpoint. Enter the new setpoint by using the number keys and accept with the ENTER key. The temperature controller will begin using the new setpoint immediately upon the press of the ENTER key.

**Note:** The sign of temperature is changed by using the arrow keys when the cursor is under the  $\pm$  character. This is possible only with temperature bath controllers.

## V.5 Data


The Review Data menu allows data review of the most recently completed speed set program. Data may be reviewed on the DV-III+ screen, on a serial printer, or both. The procedure for reviewing data after a test is as follows:

- 1) Run a DV-III or B.E.V.I.S. program.
- 2) Press the OPTION/TAB key.
- 3) Press the NUMBER 5 key to display the Review Data screen.

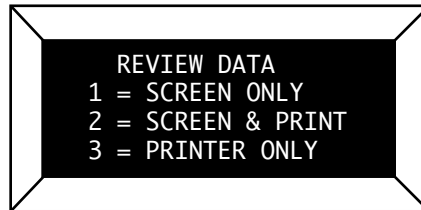

**Figure 62**

- 4) Select one of the review modes: press “NUMBER 1” for DV-III+ screen only, “NUMBER 2” for DV-III+ screen and printer, or “NUMBER 3” for printer only. Note that you should select options 2 and 3 only if the printer is connected and “on-line.”
- 5) If NUMBER 1 or NUMBER 2 were pressed, data from the first step is either displayed on the screen or displayed and printed. Press the ENTER key to scroll through and display/print the remaining step data.

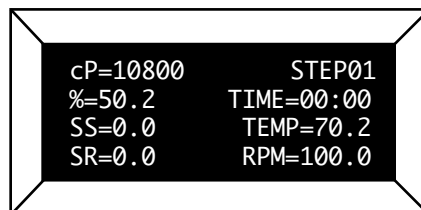

**Figure 63**

- 6) If the NUMBER 3 key (printer only) is selected, the DV-III+ prompts to “**READY PRINTER PRESS ANY KEY**” or “**PRINTER IS NOT READY TURN ON/PRESS ONLINE**”. A sample of printed output is shown in Figure 64.

○	Sample Name: _____	○
○	Operator Name: _____	○
○	Date: 06/21/91 Time: 10:06 Model: 2R Spindle: 31	○
○	#01 RPM=30.0 %=18.2 cP=970 D/CM2=100 1/SEC=10.2 T=22.1C TIME=00:05	○
○	#02 RPM=40.0 %=24.3 cP=972 D/CM2=132 1/SEC=13.6 T=22.1C TIME=00:05	○
○	#03 RPM=50.0 %=30.3 cP=969 D/CM2=164 1/SEC=17.0 T=22.1C TIME=00:05	○
○	#04 RPM=60.0 %=36.3 cP=967 D/CM2=197 1/SEC=20.4 T=22.1C TIME=00:05	○
○	#05 RPM=70.0 %=42.3 cP=966 D/CM2=230 1/SEC=23.8 T=22.1C TIME=00:05	○

**Figure 64**

**Review data notes:**

- 1) The DV-III+ stores data from the most recent speed set program test. When a subsequent speed set program is run, any test data in memory will be overwritten and the previous data will be lost.
- 2) If the “Printer Only” option is selected, the date, time of day, model and spindle number are printed with the test results. The DV-III+ does not store time and date when it is turned off; therefore, the time and date must be entered when printing the first time after start-up. See Section V.2.
- 3) The DV-III+ may be set-up to “handshake” (using **XOn/XOff** protocol) or not handshake when connected to a serial printer. If you opt to use the handshake mode, your printer is “on line”, all cable connections are correct and your printer is set for handshake mode, then data should appear on your printer immediately. If it does not appear, and instead you see the message: “**PRINTER IS NOT READY TURN ON/PRESS ONLINE**” then you should turn the handshake option off. See Section V.2.

**V.6 TIMED DATA COLLECTION**

The DV-III+ offers three methods of time control that may be used independently of control programs. These techniques result in a single data point collected at the end of the test. If multiple points are required, the rheometer should be connected to a printer and set up with an appropriate print interval.

◆ **TIME TO TORQUE**

The DV-III+ will record the amount of time required to reach the specified % torque value at a single speed. Torque and speed are input using the number keys and OPTIONS/TAB key. The test will begin immediately upon the press of the ENTER key.

The test will end when the specified torque level is reached (in either an upward or downward direction). The data will be displayed as shown in Figure 65. The SELECT DISPLAY key can be used to view all measurement parameters. The PRINT key can be used to send a single data line to the printer.

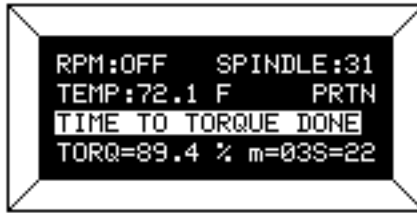


Figure 65

Note: If continuous printing was used during the test, it will be suspended when viewing the test data or entering time to torque parameters.

#### ◆ TIMED STOP

The DV-III+ will operate at a single RPM for a specified period of time. Time and speed are entered using the number keys and OPTIONS/TAB key. The test will begin immediately upon the press of the ENTER key.

The test will end when the specified time interval has elapsed. The data will be shown as displayed in Figure 66. The SELECT DISPLAY key can be used to view all measurement parameters. The PRINT key can be used to send a single dataline to the printer.

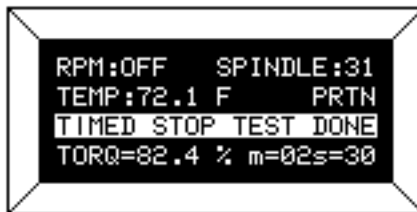


Figure 66

Note: Continuous print mode may be used in conjunction with Timed Stop to print data throughout the time period. Printing will occur only during the test.

#### ◆ TIMED AVERAGE

The DV-III+ will collect a specified number of data points over a period of time and present the average reading (arithmetic average). Time, speed and the number of data points are entered using the number keys and OPTIONS/TAB key. The test will begin immediately upon the press of the ENTER key. **Non-averaged data will be displayed during the test.**

The test will end when the specified time interval has elapsed. The averaged data point will be shown as displayed in Figure 67. The SELECT DISP key can be used to view the average value of all measurement parameters. The PRINT key can be used to send a single averaged dataline to the printer.



DATA AVERAGING  
300 READINGS WERE  
AVG'D OVER 02:30 M:S  
AVERAGE TORQUE=90.3

**Figure 67**

**Note:** Printing during a Timed Average test will show non-averaged data. Averaged data is available only at the conclusion of the test. Printing the averaged data must be done with the PRINT key at the end of the test.

## VI. Rheoloader User's Manual

B.E.V.I.S. (Brookfield Engineering Rheometer Instruction Set) is a scripting language developed at Brookfield Engineering Laboratories that allows for the creation of flexible programs to control our line of Rheometers. In the case of the DV-III+ Rheometer, programs are created then loaded into the Rheometer using the RheoLoader software.

Some features of the scripting language are:

- Repeatedly run the same program for quality control purposes.
- Wait for a prevailing condition before continuing with the program (i.e. torque value, a temperature value, a key press, etc.).
- Run the Rheometer at any of the speeds in the Custom Speed menu.
- Display messages to the screen or an attached printer to aid in operator usability.
- An internal clock that keeps time between each printed data line (this time is displayed as the last parameter on each printed line) providing a consistent time base for the collected data.

### VI.1 Description of B.E.V.I.S. Commands

Command Code	Required Parameter	Command Description
WTI	Time (MM:SS)	The program waits at this step until the specified time elapses.
WPT	% Torque value (%)	The program waits at this step until the current % torque equals the specified value.
WTP	Temperature value (°C)	The program waits at this step until the current temperature equals the specified value.
WKY	16 character (or less) text message	The specified message is displayed on the top line of the DV-III+ display while <b>PRESS A KEY</b> is displayed on the bottom line of the DV-III+. The program waits at this step until a viscometer key is pressed. While waiting at this step, the viscometer produces a beep every few seconds to remind the operator that a keypress is required to continue. If a print interval was enabled (see SPI) at the time this command is executed, the data print timer continues to count up. If the print interval elapses and a key has not yet been pressed, a line of data displaying the time since the last data print is printed as soon as a key is pressed.
SSN	Speed value (RPM)	The DV-III+ begins rotating at the specified speed. This can be any of the speeds listed in the Speed list of the DVLoader software. These speeds are the same as those listed in the Custom Speeds list in the viscometer's Options menu.
SPI	Time (MM:SS)	The DV-III+ begins printing data to the selected printer (serial or parallel; as selected in the DV-III+ menus) at the rate specified. MM:SS is minutes:seconds.
SSP	Two digit spindle code	Calculations of viscosity, shear stress, and shear rate are performed based on the specified spindle code. This command overrides the spindle currently entered via the keypad on the DV-III+.
STM	Temperature value	Set and control to the specified temperature if a Brookfield Engineering Labs. temperature controller is attached to the rheometer.
STZ	N/A	Sets the data print timer clock back to zero.
PDN	N/A	The DV-III+ immediately prints a data string to the selected printer (serial or parallel; as selected in the DV-III+ menus).
PLN	16 character (or less) text message	The DV-III+ prints the specified message to the selected printer (serial or parallel; as selected in the DV-III+ menus).

**NOTE: MM:SS is minutes:seconds**

By using various combinations of the above commands, programs are created that automatically control the viscometer and collect data (via an attached printer) from the DV-III+ Rheometer.

1	SSN	50.0
2	WTI	00:30
3	SPI	00:05
4	WPT	75 .0
5	SSN	25.0
6	WTI	00:10
7	WPT	75 .0
8	SSN	0.0
9		
10		
11		
12		
13		
14		
15		
16		
17		
18		
19		
20		
21		
22		
23		
24		
25		

The RheoLoader software is a Windows 95 (or above) based program used to create, save, print and download B.E.V.I.S. programs to the DV-III+ Rheometer. Start the software by clicking on its associated icon or by clicking the Start button; select Run; enter the name of the program to execute {Rheoad.exe}; then click OK.

This is the grid where the operator programs are created. It is used to view and edit the B.E.V.I.S. programs. When the software starts, an empty grid appears on the left of the screen. Highlight a command in the list box to the right of this grid, then click on the Insert button to insert the command into the highlighted line of the grid. This same insertion task can also be accomplished by double-clicking on the appropriate command in the list box to the right.


Before being permitted to insert another command, the parameter for the previous command in the grid must be entered if one is required. This parameter is entered into the last column of the grid.



Select the COM (RS-232) port the Programmable DV-III+ Rheometer is connected to from the COM Port drop down list.



Click the Open File button to load existing B.E.V.I.S. programs.

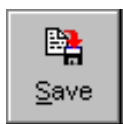

Click the Save File button to save the B.E.V.I.S. program displayed in the grid.



Click the Print button to print the B.E.V.I.S. program displayed in the grid.



Click the Insert button to insert the B.E.V.I.S. command selected in the Commands box into the selected row in the program grid.



Click the Delete button to delete the command in the selected row of the program grid.



Click the Up button to move the command in the selected row of the program grid up one row.



Click the Down button to move the command in the selected row of the program grid down one row.



Click the Clear button to clear the grid of all commands. Once cleared, the commands cannot be retrieved.

	<b>WTI</b>	<b>Wait for time interval</b>
	<b>WPT</b>	<b>Wait for % torque</b>
	<b>WTP</b>	<b>Wait for temperature</b>
	<b>WKY</b>	<b>Wait for a keypress</b>
	<b>SSN</b>	<b>Set viscometer speed</b>
	<b>SSP</b>	<b>Set spindle</b>
	<b>STZ</b>	<b>Reset data timer</b>
	<b>SPI</b>	<b>Set print interval</b>
	<b>PDN</b>	<b>Print data point now</b>
	<b>PLN</b>	<b>Print text now</b>

This list box displayed the commands available for creating programs. As previously stated, clicking on the Insert button inserts the highlighted command (WTI in this case) into the selected line in the program grid. Double-clicking on a line in this list box also inserts the command into the grid.

The icons to the left of the command descriptions indicate the type of command:



A command to wait for a condition.



A command to set a program parameter.



A command to send information to an attached printer.

Before downloading a program to the Rheometer, ensure the following have been checked:

- The appropriate cable (BEL Part No. DVP-80) is connected between the selected COM port of the PC and the Rheometer.
- The DV-III+ is at the download screen: OPTIONS/DOWNLOAD A PROGRAM/LOAD TO SLOT#x where x is slot 1,2,3, or 4.
- With the LOAD TO SLOT #x screen displayed, choose a store slot using the DV-III+ arrow keys, then press the ENTER/AUTORANGE key on the Rheometer. If after five seconds, the Rheometer cannot communicate with the RheoLoader program, the B.E.V.I.S. ERROR NO PC ATTACHED message is displayed. If a connection is established, the Download button on the PC software becomes enabled, and the DV-III+ screen displays DOWNLOADING

B.E.V.I.S. PROGRAM TO SLOT #1. Click on this button to download the displayed program to the DV-III+. When the download is complete, the DV-III+ displays DOWNLOAD DONE TO EXIT PRESS A KEY.

At this point, the program in the DV-III+ can be printed and/or run from the Rheometer.

Click on this button to exit the RheoLoader software.

## VI. 2 Example Programs

The following example programs can also be found on the RheoLoader disk that was included with the DV-III+ Rheometer:

Program 1: Pre-shear

Command	Command Description	Parameter	Comments
PLN	Print text now	Preshearing now	print user message
SSN	Set viscometer speed	50.0	run at 50 RPM
WPT	Wait for % torque	90.0	wait until 90% torque is reached
PLN	Print text now	Collecting data	print user message
SPI	Set print interval	00:10	begin printing data at 10 second intervals
SSN	Set viscometer speed	10.0	run at 10 RPM
WTI	Wait for time interval	01:40	wait at this step for 1 minute and 40 seconds, effectively printing 10 data lines

Program 2: For use with an external temperature controller

Command	Command Description	Parameter	Comments
STM	Set temperature	40.0	Set control value to 40°C
WTP	Wait for temperature	40.0	wait until temperature = 40°C (as an example, a Thermosel/Controller can be used for temperature control)
WTI	Wait for time interval	05:00	soak time; allow temperature to settle
SSN	Set viscometer speed	25.0	run at 25 RPM
SPI	Set print interval	00:30	begin printing data at 30 second intervals
WTI	Wait for time interval	06:00	wait at this step for 6 minutes, effectively printing 12 data lines

Program 3: Wait for cure

Command	Command Description	Parameter	Comments
SSP	Set spindle	31	set to a number 31 spindle
SSN	Set viscomter speed	100.0	run at 100 RPM
SPI	Set print interval	00:05	begin printing data at 5 second intervals
WPT	Wait for % torque	85.0	wait until % torque = 85; a curing cycle

Program 4: Spring relax

Command	Command Description	Parameter	Comments
WKY	Wait for a key press	Wind to 100%	tell operator to wind spindle until 100% torque is reached
WPT	Wait for % torque	100.0	wait until 100% torque is reached
WKY	Wait for a key press	Release spindle	tell operator to release the spindle
SPI	Set print interval	00:01	begin printing data at 1 second intervals
WPT	Wait for % torque	0.0	wait for spindle to completely unwind to 0% torque

Program 5: Variable speed

Command	Command Description	Parameter	Comments
SSN	Set viscometer speed	5.0	run at 50 RPM
WTI	Wait for time interval	00:10	wait for 10 seconds
PDN	Print data point now		print one data point
SSN	Set viscometer speed	10	run at 10 RPM
WTI	Wait for time interval	00:10	wait for 10 seconds
PDN	Print data point now		print one data point
SSN	Set viscometer speed	20	run at 20 RPM
WTI	Wait for time interval	00:10	wait for 10 seconds
PDN	Print data point now		print one data point
SSN	Set viscometer speed	5.0	run at 50 RPM
WTI	Wait for time interval	00:10	wait for 10 seconds
PDN	Print data point now		print one data point

## APPENDIX A - Cone/Plate Rheometer Set-Up

This Cone/Plate version of the DV-III+ uses the same operating instruction procedures as described in this manual. However, the “gap” between the cone and the plate must be verified/adjusted before measurements are made. This is done by moving the plate (built into the sample cup) up towards the cone until the pin in the center of the cone touches the surface of the plate, and then by separating (lowering) the plate 0.0005 inch (0.013mm).

Programmable DV-III+ Cone/Plate Viscometers, S/N 50969 and higher, have an Electronic Gap Setting feature. This feature enables the user to easily find the 0.0005 inch gap setting that was established at Brookfield prior to shipment.

The following information explains how to set the Electronic Gap and verify calibration of the DV-III+ Viscometer.

### A.1 ELECTRONIC GAP SETTING FEATURES


**TOGGLE SWITCH** allows you to enable/disable the Electronic Gap Setting Feature: left position is OFF (disabled), right position is ON (enabled).

**PILOT LIGHT** is the red (LED) light; when illuminated, it means the Electronic Setting Function is sensing (enabled).

**CONTACT LIGHT** is the yellow (LED) light; when it first turns on, the “hit point” has been found.

**SLIDING REFERENCE MARKER** is used after finding the “hit point;” it is the reference for establishing the 0.0005 inch gap.




**MICROMETER ADJUSTMENT RING** is used to move the cup up or down in relation to the cone spindle. Turning the ring left (clockwise) lowers the cup; turning it right (counterclockwise) raises the cup. Each line on the ring represents one scale division and is equivalent to 0.0005 inch movement of the plate relative to the cone.

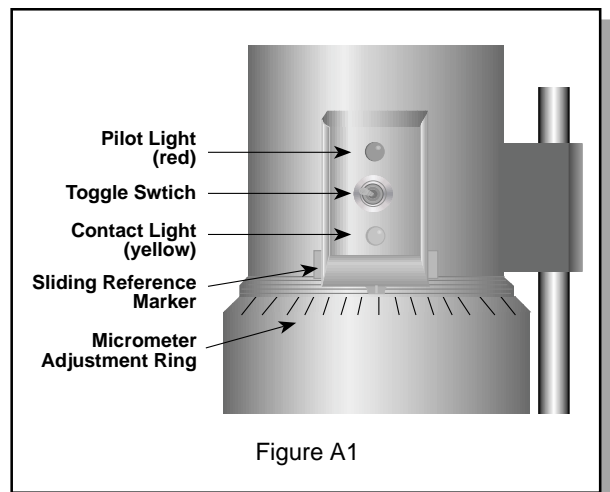

## A.2 SETUP

1. Be sure that the Viscometer is securely mounted to the Laboratory Stand, leveled and zeroed with no cone or cup attached and 0% torque is displayed.
2. **Figure A2** shows a typical water bath setup. Connect the sample cup inlet/outlet ports to the water bath inlet and outlet and set the bath to the desired test temperature. Allow sufficient time for the bath to reach the test temperature.
3. The Viscometer has been supplied with a special cone spindle(s) which contains the Electronic Gap Setting feature. The “CPE” part number designation on the cone verifies the Electronic Gap Setting feature. **Note:** The “CPE” cone or cup cannot be used with earlier DV-III+ cone/plate Viscometers (below S/N50969) which do not have the electronic gap setting feature.
4. With the motor off, thread the cone spindle by using the spindle wrench to secure the viscometer coupling nut (see **Figure A3**); gently push up on the coupling nut and hold this securely with the wrench. Thread the cone spindle by hand. **Note:** Left Hand Threads.
5. Attach the cup, taking care not to hit the cone with the cup (**Figure A4**).

The viscosity of electrically conductive fluids may be affected if readings are taken with the Electronic Gap Setting feature “on”. Be sure to shut the feature “off” before taking readings!

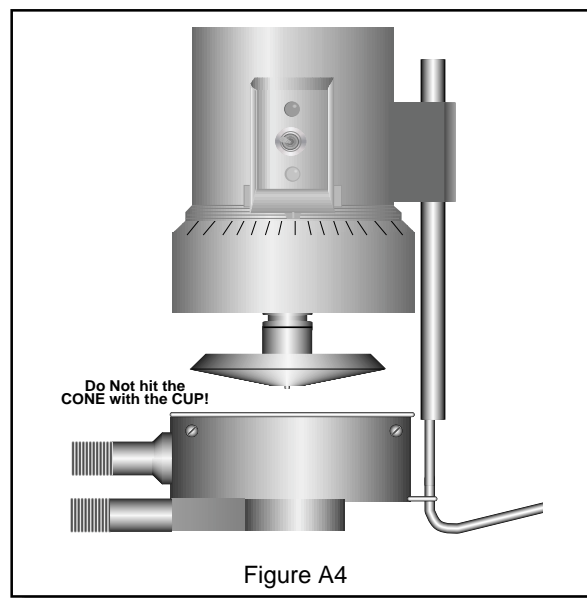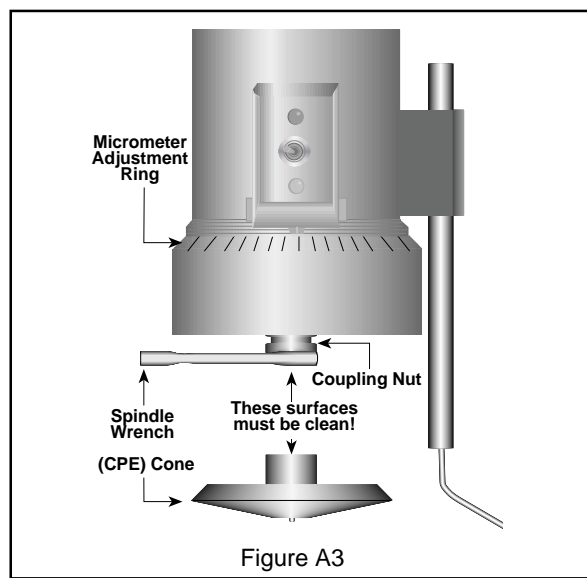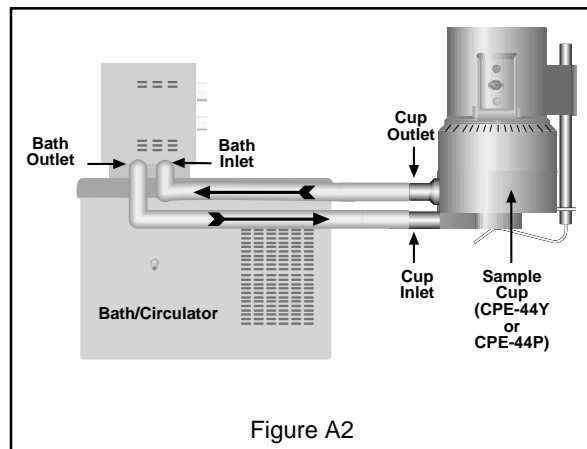

### A.3 SETTING THE GAP

1. Move the toggle switch to the right; this will turn on (enable) the Gap Setting Feature. The Pilot (red) light will be illuminated.
2. If the contact light (yellow) is illuminated, turn the micrometer adjustment ring clockwise (as you look down on the instrument) until the light is just breaking contact, i.e., flickering (see **Figure A5**).
3. If the yellow contact light is not illuminated, *slowly* turn the micrometer adjustment ring in small increments (one or two scale divisions) counter-clockwise.

Continue moving the micrometer adjustment ring *slowly* counter-clockwise until the contact light (yellow) turns on. Back off (rotate clockwise) until the light is just breaking contact, i.e., flickering.

4. Adjust the sliding reference marker, right or left, to the closest full scale division mark (see **Figure A6**).
5. Turn the micrometer adjustment ring one scale division to the left to meet the line on the sliding reference marker. **THE YELLOW CONTACT LIGHT SHOULD GO OFF.**
6. You have established the gap space needed for measurement. Now turn the toggle switch OFF (left); the red pilot light should go off.

7. Carefully remove the sample cup.

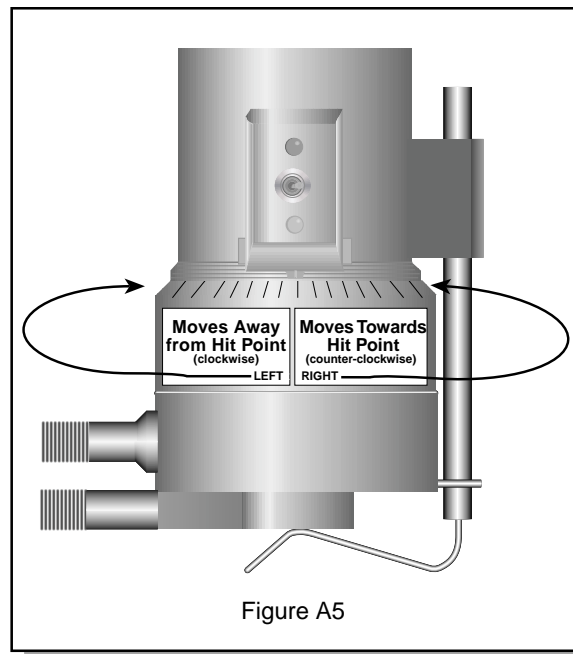


Figure A5

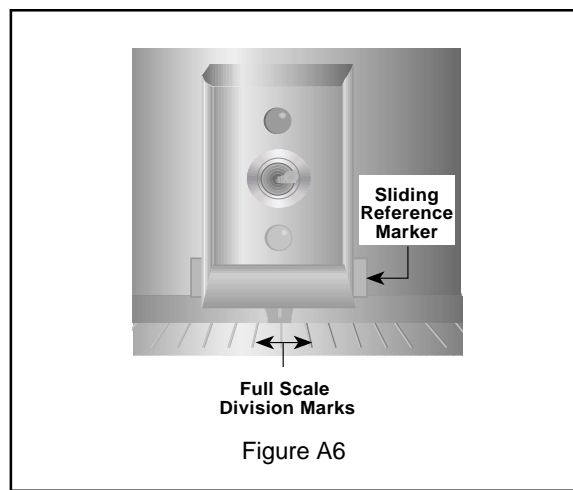


Figure A6

#### Notes

1. The cup may be removed and replaced without resetting the gap if the micrometer adjustment ring has not been moved.
2. Remove the spindle from the viscometer when cleaning.
3. Re-establish the hit point every time the spindle is attached/detached.

## A.4 VERIFYING CALIBRATION

1. Determine the appropriate sample volume. Refer to Table A1 to determine the correct sample volume required for the spindle to be utilized.
2. Select a Brookfield Viscosity Standard fluid that will give viscosity readings between 10% and 100% of full scale range. Refer to **Appendix B** for viscosity ranges of cone spindles; ranges listed apply to CPE cones.

Do not use a silicone viscosity standard fluid with a viscosity value greater than 5000 cP with a Cone/Plate. Brookfield offers a complete range of mineral oil viscosity standards suitable for use with Cone/Plates for viscosities above 5,000 cP or shear rates above 500 sec<sup>-1</sup>; see Table E2 in Appendix E for a list of available fluids.

It is best to use a viscosity standard fluid that will be close to the maximum viscosity for a given cone spindle/speed combination.

**Example:** LV DV-III+ Viscometer, Cone Spindle CPE-42, Brookfield Silicone Viscosity Standard having a viscosity of 9.7 cP at 25°C

At 60 RPM, the full scale viscosity range is 10.0 cP. Thus, the Viscometer reading should be 97% torque and 9.7 cP viscosity  $\pm$  0.197 cP. The allowable error ( $\pm$ 0.197 cP) is a combination of Viscometer accuracy and fluid tolerance (refer to **Interpretation of Calibration Test Results** in Appendix E).

3. Set the gap as described in Section A.3.
4. With the motor off, remove the sample cup and place the viscosity standard fluid into the cup.

Table A6	
Cone Part No.	Sample Volume
CPE-40	0.5 ml
CPE-41	2.0 ml
CPE-42	1.0 ml
CPE-51	0.5 ml
CPE-52	0.5 ml

Table A1

5. Attach the sample cup to the Viscometer and allow sufficient time for the sample, cup and cone to reach temperature equilibrium.
6. Turn the motor on. Set the desired speed(s). Measure the viscosity and record the reading in both % torque and centipoise (cP).

**Note:** The cone spindle must rotate at least five (5) times before a viscosity reading is taken.

7. Verify that the viscosity reading is within the allowable 1% deviation, as explained earlier, for the specific viscosity standard fluid(s) that you are using.

\* The CPE designation on the cone spindle indicates use with Electronic Gap Setting Cone/Plate Viscometers/Rheometers **only**.

## APPENDIX B - Viscosity Ranges

The table below (Universal Spindle Ranges) lists the Spindle Range Coefficients for all spindles used on DV-III+ Rheometers. Dividing the coefficient number by any of the 2,500 Rheometer speeds will give the full scale viscosity range for a Rheometer/spindle/speed combination.

### Universal Spindle Ranges

Spindle	Entry Code	SPINDLE RANGE COEFFICIENT			
		Rheometer Series			
		LV	RV	HA	HB
RV1	01	937	10,000	20,000	80,000
RV2	02	3,750	40,000	80,000	320,000
RV3	03	9,375	100,000	200,000	800,000
RV4	04	18,750	200,000	400,000	1,600,000
RV5	05	37,500	400,000	800,000	3,200,000
RV6	06	93,750	1,000,000	2,000,000	8,000,000
RV7	07	375,000	4,000,000	8,000,000	32,000,000
HA1	01	937	10,000	20,000	80,000
HA2	02	3,750	40,000	80,000	320,000
HA3	03	9,375	100,000	200,000	800,000
HA4	04	18,750	200,000	400,000	1,600,000
HA5	05	37,500	400,000	800,000	3,200,000
HA6	06	93,750	1,000,000	2,000,000	8,000,000
HA7	07	375,000	4,000,000	8,000,000	32,000,000
HB1	01	937	10,000	20,000	80,000
HB2	02	3,750	40,000	80,000	320,000
HB3	03	9,375	100,000	200,000	800,000
HB4	04	18,750	200,000	400,000	1,600,000
HB5	05	37,500	400,000	800,000	3,200,000
HB6	06	93,750	1,000,000	2,000,000	8,000,000
HB7	07	375,000	4,000,000	8,000,000	32,000,000
LV1	61	6,000	64,000	128,000	512,000
LV2	62	30,000	320,000	640,000	2,560,000
LV3	63	120,000	1,280,000	2,560,000	10,240,000
LV4	64	600,000	6,400,000	12,800,000	51,200,000
LV5	65	1,200,000	12,800,000	25,600,000	102,400,000
T-A	91	18,750	200,000	400,000	1,600,000
T-B	92	37,440	400,000	800,000	3,200,000
T-C	93	93,600	1,000,000	2,000,000	8,000,000
T-D	94	187,200	2,000,000	4,000,000	16,000,000
T-E	95	468,000	5,000,000	10,000,000	40,000,000
T-F	96	936,000	10,000,000	20,000,000	80,000,000
Spiral	70	98,400	1,050,000	2,100,000	8,400,000
ULA	00	600	6,400	12,800	51,200

## Universal Spindle Ranges

Spindle	Entry Code	SPINDLE RANGE COEFFICIENT			
		Rheometer Series			
		LV	RV	HA	HB
DIN-81	81	3,420	36,500	73,000	292,000
DIN-82	82	3,420	36,500	73,000	292,000
DIN-83	83	11,340	121,300	242,600	970,400
DIN-85	85	1,144	12,200	24,400	97,600
DIN-86	86	3,420	36,500	73,000	292,000
SC4-14	14	117,200	1,250,000	2,500,000	10,000,000
SC4-15	15	46,880	500,000	1,000,000	4,000,000
SC4-16	16	120,000	1,280,000	2,560,000	10,240,000
SC4-18	18	3,000	32,000	64,000	256,000
SC4-21	21	4,688	50,000	100,000	400,000
SC4-25	25	480,000	5,120,000	10,240,000	40,960,000
SC4-27	27	23,4000	250,000	500,000	2,000,000
SC4-28	28	46,8880	500,000	1,000,000	4,000,000
SC4-29	29	93,750	1,000,000	2,000,000	8,000,000
SC4-31	31	30,000	320,000	640,000	2,560,000
SC4-34	34	60,000	640,000	1,280,000	5,120,000
SC4-37	37	23,4000	250,000	500,000	2,000,000
CP-40	40	307	3,270	6,540	26,160
CP-41	41	1,151	12,280	24,560	98,240
CP-42	42	600	6,400	12,800	51,200
CP-51	51	4,800	51,200	102,400	409,600
CP-52	52	9,216	98,300	196,600	786,400

**Notes:** RV spindle ranges are calculated with the RV Guardleg in use.  
LV spindle ranges are calculated with the LV Guardleg in use.

**Example 1:** Determine the full scale viscosity range of the LV3 spindle running on an RV series Rheometer at 45 RPM.

Spindle Speed = 45 RPM

LV Spindle Range Coefficient for RV series Rheometer = 1,280,000

Full Scale Viscosity Range =  $\frac{1,280,000}{45} = 28,444$  cP (mPa\*s)

**Example 2:** Determine the full scale viscosity range of the LV3 spindle running on a 2xHA Rheometer at 103.5 RPM.

Spindle Speed = 103.5 RPM

LV3 Spindle Range Coefficient for 2xHA Rheometer = 2 x 2,560,000

Full Scale Viscosity Range =  $\frac{5,120,000}{103.5} = 49,468$  cP (mPa\*s)

**Note:** The maximum viscosity that should be taken from the DV-III+ Rheometer is at 100% of any full scale spindle/speed range.

## LV(#1-4) and RV,HA,HB(#1-7) Rheometers

Viscosity Range (cP)		
Viscometer	Minimum	Maximum
LVDV-III+	15	6,000,000
RVDV-III+	100	40,000,000
HADV-III+	200	80,000,000
HBDV-III+	800	320,000,000

## Small Sample Adapter

Spindle	Shear Rate (sec <sup>-1</sup> )	Viscosity (cP)			
		LVDV-III+cP	RVDV-III+cP	HADV-III+cP	HBDV-III+cP
14	0 - 80	47.0 - 1,171,000	500.0 - 12,500,000	1,000.0 - 25,000,000	4,000 - 100,000,000
15	0 - 96	19.0 - 468,650	200.0 - 5,000,000	400.0 - 10,000,000	1,600 - 40,000,000
16	0 - 58	48.0 - 1,199,700	512.0 - 12,800,000	1,024.0 - 25,600,000	4,100 - 102,400,000
18	0 - 264	1.3 - 30,000	12.8 - 320,000	25.6 - 640,000	103 - 2,560,000
21	0 - 186	1.9 - 46,865	20.0 - 500,000	40.0 - 1,000,000	160 - 4,000,000
25	0 - 44	192.0 - 4,790,000	2,050.0 - 51,200,000	4,100.0 - 102,400,000	16,400 - 409,600,000
27	0 - 68	9.4 - 234,325	100.0 - 2,500,000	200.0 - 5,000,000	800 - 20,000,000
28	0 - 56	18.8 - 468,650	200.0 - 5,000,000	400.0 - 10,000,000	1,600 - 40,000,000
29	0 - 50	37.5 - 937,300	400.0 - 10,000,000	800.0 - 20,000,000	3,200 - 80,000,000
31	0 - 68	12.0 - 300,000	128.0 - 3,200,000	256.0 - 6,400,000	1,024 - 25,600,000
34	0 - 56	24.0 - 600,000	256.0 - 6,400,000	512.0 - 12,800,000	2,048 - 51,200,000
37	0 - 72	9.4 - 234,325	100.0 - 2,500,000	200.0 - 5,000,000	800 - 20,000,000
82	0 - 258	1.4 - 10,000	15.0 - 10,000	30.0 - 10,000	120 - 10,000
83	0 - 258	4.6 - 50,000	50.0 - 50,000	100.0 - 50,000	400 - 50,000

## Thermosel

Spindle	Shear Rate (sec <sup>-1</sup> )	Viscosity (cP)			
		LVDV-III+cP	RVDV-III+cP	HADV-III+cP	HBDV-III+cP
18	0 - 264	1.3 - 30,000	12.8 - 320,000	25.6 - 640,000	103 - 2,560,000
21	0 - 186	1.9 - 46,865	20.0 - 500,000	40.0 - 1,000,000	160 - 4,000,000
27	0 - 68	9.4 - 234,325	100.0 - 2,500,000	200.0 - 5,000,000	800 - 20,000,000
28	0 - 56	18.8 - 468,650	200.0 - 5,000,000	400.0 - 10,000,000	1,600 - 40,000,000
29	0 - 50	37.5 - 937,300	400.0 - 10,000,000	800.0 - 20,000,000	3,200 - 80,000,000
31	0 - 68	12.0 - 300,000	128.0 - 3,200,000	256.0 - 6,400,000	1,024 - 25,600,000
34	0 - 56	24.0 - 600,000	256.0 - 6,400,000	512.0 - 12,800,000	2,048 - 51,200,000
81	0 - 258	1.4 - 10,000	15.0 - 10,000	30.0 - 10,000	120 - 10,000

## UL Adapter

UL Spindle	Shear Rate (sec <sup>-1</sup> )	Viscosity (cP)			
		LVDV-III+	RVDV-III+	HADV-III+	HBDV-III+
YULA-15 or 15Z	0 - 122.4	1.0 - 2,000	3.2 - 2,000	6.4 - 2,000	25.6 - 2,000

**Note:** LV DV-III+ - Minimum viscosity, 1.0 cP at 60RPM, 10% of full scale range.

**Note:** RV,HA,HB-III - Minimum viscosity; 3 cP, 6 cP and 24 cP, at 230.9 RPM, 10% of full scale range.

## Din Adapter

DAA Spindle	Shear Rate (sec-1)	Viscosity (cP)			
		LVDV-III+	RVDV-III+	HADV-III+	HBDV-III+
85	0 - 258	0.6 - 5,000	6.1 - 5,000	12.2 - 5,000	48.8 - 5,000
86	0 - 258	1.8 - 10,000	18.2 - 10,000	36.5 - 10,000	146 - 10,000
87	0 - 258	5.7 - 50,000	61.0 - 50,000	121 - 50,000	485 - 50,000

## Spiral Adapter

Spiral Spindle	Shear Rate (sec-1)	Viscosity (cP)			
		LVDV-III+	RVDV-III+	HADV-III+	HBDV-III+
SA-70	.677 - 67.7 (1 - 100 RPM)	98.4 - 98,416	1,050 - 1,050,000	2,100 - 2,100,000	8,400 - 8,400,000

## Cone/Plate Rheometer

Cone Spindle	Shear Rate (sec-1)	Viscosity (cP)			
		LVDV-III+C/P	RVDV-III+C/P	HADV-II+IC/P	HBDV-II+IC/P
CPE-40	0 - 1500	0.15 - 3,065	1.7 - 32,700	3.3 - 65,400	13.1 - 261,000
CPE-41	0 - 400	0.58 - 11,510	6.2 - 122,800	12.3 - 245,600	49.1 - 982,400
CPE-42	0 - 760	0.3 - 6,000	3.2 - 64,000	6.4 - 128,000	25.6 - 512,000
CPE-51	0 - 768	2.4 - 47,990	25.6 - 512,000	51.2 - 1,024,000	205.0 - 4,096,000
CPE-52	0 - 400	4.6 - 92,130	49.2 - 983,000	98.3 - 1,966,000	393.0 - 7,864,000

## Helipath with T-Bar Spindles

T-Bar Spindle	Viscosity (cP)			
	LVDV-III+	RVDV-III+	HADV-III+	HBDV-III+
T-A	156 - 187,460	2,000 - 200,000	4,000 - 4,000,000	16,000 - 16,000,000
T-B	312 - 374,920	4,000 - 400,000	8,000 - 8,000,000	32,000 - 32,000,000
T-C	780 - 937,300	10,000 - 10,000,000	20,000 - 20,000,000	80,000 - 80,000,000
T-D	1,560 - 1,874,600	20,000 - 20,000,000	40,000 - 40,000,000	160,000 - 160,000,000
T-E	3,900 - 4,686,500	50,000 - 50,000,000	100,000 - 100,000,000	400,000 - 400,000,000
T-F	7,800 - 9,373,000	100,000 - 100,000,000	200,000 - 200,000,000	800,000 - 800,000,000

In taking viscosity measurements with the DV-III+ Rheometer there are two considerations which pertain to the low viscosity limit of effective measurement.

- 1) Viscosity measurements should be accepted within the equivalent % Torque Range from 10% to 100% for any spindle/speed combination.
- 2) Viscosity measurements should be taken under laminar flow conditions, not under turbulent flow conditions.

The first consideration has to do with the precision of the instrument. All DV-III+ Rheometers have a

full scale range precision of (+/-) 1% of any spindle/speed combination. We discourage taking readings below 10% of range because the potential viscosity error of (+/-) 1% is a relatively high number compared to the instrument reading.

The second consideration involves the mechanics of fluid flow. All rheological measurements of fluid flow properties should be made under laminar flow conditions. Laminar flow is flow wherein all particle movement is in layers directed by the shearing force. For rotational systems, this means all fluid movement must be circumferential. When the inertial forces on the fluid become too great, the fluid can break into turbulent flow wherein the movement of fluid particles becomes random, and the flow can not be analyzed with standard math models. This turbulence creates a falsely high Rheometer reading, with the degree of non-linear increase in reading being directly related to the degree of turbulence in the fluid.

For the following geometries, we have found that an approximate transition point to turbulent flow occurs:

- 1) No. 1 LV Spindle: 15 cP at 60 RPM
- 2) No. 1 RV Spindle: 100 cP at 50 RPM
- 3) UL Adapter: 0.85 cP at around 70 RPM
- 4) SC4-18/13R: 1.25 cP at around 240 RPM

Turbulent conditions will exist in these situations whenever the RPM/cP ratio exceeds the values listed above.

## APPENDIX C - Variables in Viscosity Measurements

As with any instrument, there are variables that can affect a viscosity measurement. These variables may be related to the instrument (Rheometer) or the test fluid. Variables related to the test fluid deal with the rheological properties of the fluid, while instrument variables would include the Rheometer design and the spindle geometry system utilized.

### RHEOLOGICAL PROPERTIES

Fluids have different rheological characteristics that can be described by Rheometer measurements. We can then work with these fluids to suit our lab or process conditions.

There are two categories of fluids:

- Newtonian** - These fluids have the same viscosity at different Shear Rates (different RPM's) and are called Newtonian over the Shear Rate range they are measured.
- Non-Newtonian** - These fluids have different viscosities at different shear rates (different RPM's). They fall into two groups:
  - 1) Time Independent non-Newtonian
  - 2) Time Dependent non-Newtonian

The time dependency is the time they are held at a given Shear Rate (RPM). They are non-Newtonian, and when you change the Rheometer spindle speed, you get a different viscosity.

#### Time Independent

- Pseudoplastic** - A pseudoplastic material displays a decrease in viscosity with an increase in shear rate, and is also known as “shear thinning”. If you take Rheometer readings from a low to a high RPM and then back to the low RPM, and the readings fall upon themselves, the material is time independent pseudoplastic and shear thinning.

#### Time Dependent

- Thixotropic** - A thixotropic material has decreasing viscosity under constant shear rate. If you set a Rheometer at a constant speed, recording viscosity (cP) values over time and find that the viscosity (cP) values decrease with time, the material is thixotropic.

Brookfield publication, “*More Solutions to Sticky Problems*” includes a more detailed discussion of rheological properties and non-Newtonian behavior.

### Rheometer Related Variables

Most fluid viscosities are found to be non-Newtonian. They are Shear Rate dependent on the measurement conditions. The specifications of the Rheometer spindle and chamber geometry will affect the viscosity readings. If one reading is taken at 25 rpm, and a second at 50 rpm, the two viscosity (cP) values produced will be different because the readings were made at different shear rates. The faster the spindle speed, the higher the shear rate.

The shear rate of a given measurement is determined by: the rotational speed of the spindle, the size and shape of the spindle, the size and shape of the container used, and therefore, the distance between the container wall and the spindle surface.

A repeatable viscosity test should control or specify the following:

- 1) Test temperature
- 2) Sample container size (or spindle/chamber geometry)
- 3) Sample volume
- 4) Rheometer model
- 5) Spindle used (if using LVDV-III+ (#1-4) or RVDV-III+ (#1-7) attach the guard leg)
- 6) Test speed or speeds (or the shear rate)
- 7) Length of time or number of spindle revolutions to record viscosity.

## APPENDIX D - Spindle and Model Codes

Each spindle has a two digit code which is entered using the SPDL key on the DV-III+ key pad. The entry code allows the DV-III+ to calculate Viscosity, Shear Rate and Shear Stress values.

Each spindle has two constants which are used in these calculations. The Spindle Multiplier Constant (SMC) used for viscosity calculations, and the Shear Rate Constant (SRC), used for shear rate and shear stress calculations. Note that where SRC = 0, no shear rate/shear stress calculations are done and the data displayed is zero (0) for these functions.

*Table D1*

SPINDLE	ENTRY CODE	SMC	SRC
RV1	01	1	0
RV2	02	4	0
RV3	03	10	0
RV4	04	20	0
RV5	05	40	0
RV6	06	100	0
RV7	07	400	0
HA1	01	1	0
HA2	02	4	0
HA3	03	10	0
HA4	04	20	0
HA5	05	40	0
HA6	06	100	0
HA7	07	400	0
HB1	01	1	0
HB2	02	4	0
HB3	03	10	0
HB4	04	20	0
HB5	05	40	0
HB6	06	100	0
HB7	07	400	0
LV1	61	6.4	0
LV2	62	32	0
LV3	63	128	0
LV4	64	640	0
LV5	65	1280	0
T-A	91	20	0
T-B	92	40	0
T-C	93	100	0
T-D	94	200	0
T-E	95	500	0
T-F	96	1000	0
ULA	00	0.64	1.223
SC4-14	14	125	0.4
SC4-15	15	50	0.48
SC4-16	16	128	0.2929
SC4-18	18	3.2	1.32
SC4-21	21	5	0.93
SC4-25	25	512	0.22

*Table D1 (continued)*

SPINDLE	ENTRY CODE	SMC	SRC
SC4-27	27	25	0.34
SC4-28	28	50	0.28
SC4-29	29	100	0.25
SC4-31	31	32	0.34
SC4-34	34	64	0.28
SC4-37	37	25	0.36
CP40	40	0.327	7.5
CP41	41	1.228	2
CP42	42	0.64	3.8
CP51	51	5.12	3.84
CP52	52	9.83	2
Spiral Adapter	70	105	0.677
Thermosel DIN spindle	81	3.7	1.29
SSA DIN spindle for 13R or 13RP chamber	82	3.75	1.29
SSA DIN spindle for 7R or 7RP chamber	83	12.09	1.29
ULA DIN spindle	85	1.22	1.29
ULA DIN spindle	86	3.65	1.29
ULA DIN spindle	87	12.13	1.29

Table D2 lists the model codes and spring torque constants for each Rheometer model.

*Table D2*

MODEL	TK	MODEL CODE ON DV-III SCREEN
LVDV-III+	0.09373	LV
2.5xLVDV-III+	0.2343	4L
5xLVDV-III+	0.4686	5L
1/4 RVDV-III+	0.25	1R
1/2 RVDV-III+	0.5	2R
RVDV-III+	1	RV
HADV-III+	2	HA
2xHADV-III+	4	3A
2.5xHADV-III+	5	4A
HBDV-III+	8	HB
2xHBDV-III+	16	3B
2.5xHBDV-III+	20	4B
5xHBDV-III+	40	5B

The full scale viscosity range for any DV-III+ model and spindle may be calculated using the equation:

$$\text{where: Full Scale Viscosity Range [cP]} = \text{TK} * \text{SMC} * \frac{10,000}{\text{RPM}}$$

**TK** = DV-III+ Torque Constant from Table D2

**SMC** = Spindle Multiplier Constant listed in Table D1

The Shear Rate calculation is:

$$\text{Shear Rate (1/sec)} = \text{SRC} * \text{RPM}$$

where:

**SRC** = Shear Rate Constant from Table D1

### Using Non-standard spindles with DV-III+ and RHEOCALC Software

Spindle Entry 99 allows entry of spindle constants which the DV-III+ will use to calculate Viscosity, Shear Rate and Shear Stress for spindles in boundary conditions other than the 600ml beaker or specified chamber.

The spindles must conform to geometries that allow for mathematical calculations of Shear Rate and Shear Stress i.e. coaxial cylinder.

Two constants are required:

- a) **SMC** (Spindle Multiplier Constant) which is used to calculate cP value.
- b) **SRC** (Shear Rate Constant) which is used to calculate Shear Rate and Shear Stress.

If SRC=0 (Example, LV #1-4, RV,HA,HB #1-7 spindles), Shear Rate and Shear Stress values are not calculated and can not be displayed.

The SMC value for a spindle may be calculated as follows:

- a) For new spindle conditions you calculate the SMC using a Newtonian fluid of known viscosity (Brookfield Viscosity Standard). This is done in the container with the new dimensions at the controlled temperature specified for the viscosity standard fluid.
- b) The new full scale viscosity range is calculated for a selected RPM:

$$R1 = \frac{100 n}{Y}$$

where:

R1 = The new full scale viscosity range  
n = The viscosity (in cP) of the Newtonian fluid  
Y = The Torque % reading at the selected RPM

- c) The SMC Value is then calculated:

$$\text{SMC} = \frac{[\text{Full Scale Viscosity Range (cP)}] * [\text{Selected RPM}]}{\text{TK} * 10,000}$$

where:

**TK** = DV-III+ Torque Constant from Table D2

The SRC value may be calculated for cylindrical spindle geometry using the following equation:

$$\text{Shear Rate (1/sec)} = \frac{2\omega R_c^2 R_b^2}{\chi^2 [R_c^2 - R_b^2]}$$

where:

$R_c$  = Radius of the container (in centimeters)

$R_b$  = Radius of the spindle (in centimeters)

$\chi$  = Radius at which the shear rate is to be calculated (normally the same value as  $R_b$ ; in centimeters)

$\omega$  = Angular velocity of the spindle (Rad/Sec)

$$\omega = \frac{2\pi}{60} * N$$

N = Spindle speed in RPM

SMC and SRC values are entered in RHEOCALC software. See the HELP file for details.

## APPENDIX E - Calibration Procedures

The accuracy of the DV-III+ is verified using viscosity standard fluids which are available from Brookfield Engineering Laboratories or your local Brookfield agent. Viscosity standards are Newtonian, and therefore, have the same viscosity regardless of spindle speed (or shear rate). Viscosity standards, calibrated at 25°C, are shown in Table E1.

Container size: For Viscosity Standards < 30,000 cP, use a 600 ml Low Form Griffin Beaker having a working volume of 500 ml.

For Viscosity Standards ≥ 30,000 cP, use the fluid container.

Inside Diameter: 3.25"(8.25cm)

Height: 4.75"(12.1cm)

Note: Container may be larger, but may not be smaller.

Temperature: As stated on the fluid standard label: (+/-) 0.1°C

Conditions: The DV-III+ should be set according to the operating instructions. The water bath should be stabilized at test temperature. Rheometers with the letters "LV" or "RV" in the model designation should have the guard leg attached.

*Table E1*

<b>Normal 25°C Standard Fluids</b>		<b>High Temperature Standard Fluids</b>
<b><u>Viscosity (cP)</u></b>	<b><u>Viscosity (cP)</u></b>	<b><u>Three Viscosity/Temperatures**</u></b>
5	5,000	HT-30,000
10	12,500	HT-60,000
50	30,000	HT-100,000
100	60,000	
500	100,000	**25°C, 93.3°C, 149°C
1,000		Refer to Brookfield catalog for more information.

*Table E2*

<b>MINERAL OIL VISCOSITY STANDARD FLUIDS</b>	
<b>BEL Part No.</b>	<b>Viscosity (cP) 25°C</b>
B31	31
B210	210
B750	750
B1400	1,400
B2000	2,000
B11000	11,000
B20000	20,000
B80000	80,000
B200000	200,000
B420000	420,000

## Brookfield Viscosity Standard Fluid - General Information

We recommend that Brookfield Viscosity Standard Fluids be replaced on an annual basis, one year from date of initial use. These fluids are pure silicone and are not subject to change over time. However, exposure to outside contaminants through normal use requires replacement on an annual basis. Contamination may occur by the introduction of solvent, standard of different viscosity or other foreign material.

Viscosity Standard Fluids may be stored under normal laboratory conditions. Disposal should be in accordance with state, local and federal regulations as specified on the material safety data sheet.

Brookfield Engineering Laboratories does not recertify Viscosity Standard Fluids. We will issue duplicate copies of the Certificate of Calibration for any fluid within two years of the purchase date.

Brookfield Viscosity Standard Fluids are reusable provided they are not contaminated. Normal practice for usage in a 600 ml beaker is to return the material from the beaker back into the bottle. When using smaller volumes in accessories such as Small Sample Adapter, UL Adapter or Thermosel, the fluid is normally discarded.

### Calibration Procedure for LV(#1-4) and RV,HA,HB(#1-7) Brookfield spindles:

- 1) Place the viscosity standard fluid (in the proper container) into the water bath.
- 2) Lower the DV-III+ into measurement position (with guard leg if LV or RV series Rheometer is used).
- 3) Attach the spindle to the Rheometer. If you are using a disk shaped spindle, avoid trapping air bubbles beneath the disk by first immersing the spindle at an angle, and then connecting it to the Rheometer.
- 4) The viscosity standard fluid, together with the spindle, should be immersed in the bath for a minimum of 1 hour, stirring the fluid periodically, prior to taking measurements.
- 5) After 1 hour, check the temperature of the viscosity standard fluid with an accurate thermometer.
- 6) If the fluid is at test temperature ( $\pm 0.1^{\circ}\text{C}$  of the specified temperature, normally  $25^{\circ}\text{C}$ ), measure the viscosity and record the Rheometer reading.

**Note:** The spindle must rotate at least five (5) times before readings are taken.

- 7) The viscosity reading should equal the cP value on the viscosity fluid standard to within the combined accuracies of the Rheometer and the standard (as discussed in the section entitled, Interpretation of Calibration Test Results).

### Calibration Procedure for a Small Sample Adapter

When a Small Sample Adapter is used, the water jacket is connected to the water bath and the water is stabilized at the proper temperature:

- 1) Put the proper amount of viscosity standard fluid into the sample chamber. The amount varies with each spindle/chamber combination. (Refer to the Small Sample Adapter instruction manual.)
- 2) Place the sample chamber into the water jacket.

- 3) Put the spindle into the test fluid and attach the extension link, coupling nut and free hanging spindle (or directly attach the solid shaft spindle) to the DV-III+.
- 4) Allow 30 minutes for the viscosity standard, sample chamber and spindle to reach test temperature.
- 5) Measure the viscosity and record the Rheometer reading.

**Note:** The spindle must rotate at least five (5) times before a viscosity reading is taken.

### **Calibration Procedure for a Thermoel System**

When a Thermoel System is used, the controller stabilizes the Thermo Container at the test temperature. A two-step process is recommended for the Thermoel:

- A) Evaluate the calibration of the Viscometer alone according to the procedure outlined in Appendix E, "Calibration Procedure for LV (#1-4) and RV, HA, HB (#1-7) Brookfield spindles."
- B) Evaluate the Viscometer with Thermoel according to the procedure listed below:
  - a) Put the proper amount of HT viscosity standard fluid into the HT-2 sample chamber. The amount varies with the spindle used. (Refer to the Thermoel instruction manual).
  - b) Place the sample chamber into the Thermo Container.
  - c) Put the spindle into the test fluid and attach the extension link, coupling nut and free hanging spindle (or directly attach the solid shaft spindle) to the DV-III+.
  - d) Allow 30 minutes for the viscosity standard, sample chamber and spindle to reach test temperature.
  - e) Measure the viscosity and record the Rheometer reading.

**Note:** The spindle must rotate at least five (5) times before a viscosity reading is taken.

### **Calibration Procedure using UL or DIN UL Adapters**

When a UL or DIN UL Adapter is used, the water bath is stabilized at the proper temperature:

- 1) Put the proper amount of viscosity standard fluid into the UL Tube. (Refer to the UL Adapter instruction manual).
- 2) Attach the spindle (with extension link and coupling nut) onto the DV-III+.
- 3) Attach the tube to the mounting channel.
- 4) Lower the tube into the water bath reservoir, or if using the ULA-40Y water jacket, connect the inlet/outlets to the bath external circulating pump.
- 5) Allow 30 minutes for the viscosity standard, sample chamber and spindle to reach test temperature.
- 6) Measure the viscosity and record the Rheometer reading.

**Note:** The spindle must rotate at least five (5) times before a viscosity reading is taken.

## Calibration Procedure using a Helipath Stand and T-Bar Spindles

When a Helipath Stand and T-Bar spindles are used:

Remove the T-bar spindle and select a standard LV(#1-4) or RV,HA,HB(#1-7) spindle. Follow the procedures for LV(#1-4) and RV,HA,HB (#1-7) Brookfield spindles outlined above.

T-Bar spindles should not be used for verifying calibration of the DV-III+ Rheometer.

## Calibration Procedure for Spiral Adapter

- 1) Place the viscosity standard fluid (in the proper container) into the water bath.
- 2) Attach the spindle to the viscometer. Attach chamber (SA-1Y) and clamp to the viscometer.
- 3) Lower the DV-III+ into measurement position. Operate the viscometer at 50 or 60 RPM until the chamber is fully flooded.
- 4) The viscosity standard fluid, together with the spindle, should be immersed in the bath for a minimum of 1 hour, stirring the fluid periodically (operate at 50 or 60 RPM periodically), prior to taking measurements.
- 5) After 1 hour, check the temperature of the viscosity standard fluid with an accurate thermometer.
- 6) If the fluid is at test temperature ( $\pm 0.1^{\circ}\text{C}$  of the specified temperature, normally  $25^{\circ}\text{C}$ ), measure the viscosity and record the viscometer reading.

**Note:** The spindle must rotate at least five (5) times for one minute, whichever is greater before readings are taken.

- 7) The viscosity reading should equal the cP value on the viscosity fluid standard to within the combined accuracies of the viscometer and the standard (as discussed in the section entitled, Interpretation of Calibration Test Results).

## Calibration Procedure for Cone/Plate Viscometers:

- 1) Follow the above procedures for mechanically adjusting the setting of the cone spindle to the plate.
- 2) Refer to Appendix A, Table A1, and determine the correct sample volume required for the selected spindle.
- 3) Select a viscosity standard fluid that will give viscosity readings between 10% and 100% of full scale range. Refer to Appendix B for viscosity ranges of cone spindles. Consult with Brookfield or an authorized dealer to determine which fluid is appropriate.

It is best to use a viscosity standard fluid that will be close to the maximum viscosity for a given cone spindle/speed combination.

**Example:** LVDV-III Viscometer, Cone CP-42, Fluid 10  
Having a viscosity of 9.7 cP at  $25^{\circ}\text{C}$

At 60 RPM, the full scale viscosity range is 10.0 cP. Thus, the Rheometer reading should be 97% torque and 9.7 cP viscosity  $\pm 0.197$  cP (0.1 cP for the rheometer plus 0.97 cP for the fluid). The accuracy is a combination of Rheometer and fluid tolerance (refer to Interpretation of Calibration Test Results).

- 4) With the viscometer stopped, remove the sample cup and place the viscosity standard fluid into the cup, waiting 10 minutes for temperature equilibrium.
- 5) Connect the sample cup to the Rheometer. Allow sufficient time for temperature to reach equilibrium. Typically 15 minutes is the maximum time that you must wait. Less time is required if spindle and cup are already at test temperature.
- 6) Measure the viscosity and record the Rheometer reading in both % torque and centipoise (cP).

- Notes:**
- 1) The spindle must rotate at least five (5) times before a viscosity reading is taken.
  - 2) The use of Brookfield Viscosity Standard fluids in the range of 5 cP to 5000 cP is recommended for cone/plate instruments. Please contact Brookfield Engineering Laboratories or an authorized dealer if your calibration procedure requires more viscous standards.
  - 3) Select a viscosity standard fluid that will give viscosity readings between 10% and 100% of full scale range. Refer to Appendix B for viscosity ranges of cone spindles. Do not use a silicone viscosity standard fluid with a viscosity value greater than 5000 cP with a Cone/Plate Viscometer. Brookfield offers a complete range of mineral oil viscosity standards suitable for use with Cone/Plate Viscometers as shown in Table E2. Consult with Brookfield or an authorized dealer to determine which fluid is appropriate.

### Interpretation of Calibration Test Results:

When verifying the calibration of the DV-III+, the instrument and viscosity standard fluid error must be combined to calculate the total allowable error.

The DV-III+ is accurate to (+/-) 1% of any full scale spindle/speed viscosity range.

Brookfield Viscosity Standards Fluids are accurate to (+/-) 1% of their stated value.

**Example:** Calculate the acceptable range of viscosity using RVDV-III with RV-3 Spindle at 2 RPM; Brookfield Standard Fluid 12,500 with a viscosity of 12,257 cP at 25°C:

- 1) Calculate full scale viscosity range using the equation:

$$\text{Full Scale Viscosity Range [cP]} = \text{TK} * \text{SMC} * \frac{10,000}{\text{RPM}}$$

Where:

$$\begin{aligned} \text{TK} &= 1.0 \text{ from Table D2} \\ \text{SMC} &= 10 \text{ from Table D1} \end{aligned}$$

$$\text{Full Scale Viscosity Range} = \frac{1 * 10 * 10,000}{2} = 50,000 \text{ cP}$$


The viscosity is accurate to (+/-) 500 cP (which is 1% of 50,000)

- 2) The viscosity standard fluid is 12,257 cP. Its accuracy is (+/-)1% of 12,257 or (+/-)122.57 cP.
- 3) Total allowable error is (122.57 + 500) cP = (+/-) 622.57 cP.
- 4) Therefore, any viscosity reading between 11,634.4 and 12,879.6 cP indicates that the Rheometer is operating correctly. Any reading outside these limits may indicate a Rheometer problem. Contact the Brookfield technical sales department or your local Brookfield dealer/distributor with test results to determine the nature of the problem.

## The Brookfield Guardleg

The *guard leg* was originally designed to protect the spindle during use. The first applications of the Brookfield Viscometer included hand held operation while measuring fluids in a 55 gallon drum. It is clear that under those conditions the potential for damage to the spindle was great. Original construction included a sleeve that protected the spindle from side impact. Early RV guard legs attached to the dial housing and LV guard legs attached to the bottom of the pivot cup with a twist and lock mechanism.

The current guard leg is a band of metal in the shape of the letter U with a bracket at the top that attaches to the pivot cup of a Brookfield Viscometer/Rheometer. Because it must attach to the pivot cup, the guard leg cannot be used with a Cone/Plate instrument. A guard leg is supplied with all LV and RV series instruments, but not with the HA or HB series. It's shape (shown in *Figure 1*) is designed to accommodate the spindles of the appropriate spindle set; therefore, the RV guard leg is wider than the LV due to the large diameter of the RV #1 spindle. They are not interchangeable.

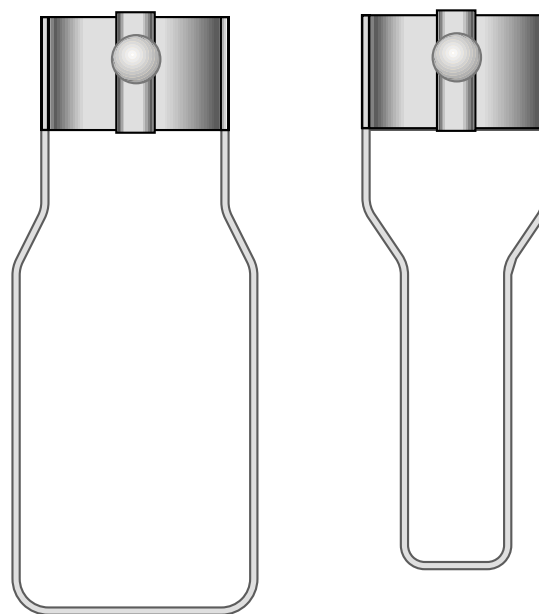

*RV Guardleg*

*LV Guardleg*

The calibration of the Brookfield Viscometer/Rheometer is determined using a 600 ml Low Form Griffin Beaker. The calibration of LV and RV series instruments includes the guard leg. The beaker wall (for HA/HB instruments) or the guard leg (for LV/RV instruments) define what is called the “outer boundary” of the measurement. The spindle factors for the LV, RV, and HA/HB spindles were developed with the above boundary conditions. The spindle *factors* are used to convert the instrument torque (expressed as the dial reading or %Torque value) into centipoise. Theoretically, if measurements are made with different boundary conditions, e.g., without the guard leg or in a container other than 600 ml beaker, then the spindle factors found on the Factor Finder cannot be used to accurately calculate an absolute viscosity. Changing the boundary conditions does not change the viscosity of the fluid, but it does change how the instrument torque is converted to centipoise. Without changing the spindle factor to suit the new boundary conditions, the calculation from instrument torque to viscosity will be incorrect.

Practically speaking, the guard leg has the greatest effect when used with the #1 & #2 spindles of the LV and RV spindle sets. Any other LV (#3 & #4) or RV (#3 - #7) spindle can be used in a 600 ml beaker with or without the guard leg to produce correct results. The HA and HB series Viscometers/Rheometers are not supplied with guard legs in order to reduce the potential problems when measuring high viscosity materials. HA/HB spindles #3 through #7 are identical to those spindle numbers in the RV spindle set. The HA/HB #1 & #2 have slightly different dimensions than the corresponding RV spindles. This dimensional difference allows the factors between the RV and HA/HB #1 & #2 spindles to follow the same ratios as the instrument torque even though the boundary conditions are different.


The recommended procedures of using a 600 ml beaker and the guard leg are difficult for some customers to follow. The guard leg is one more item to clean. In some applications the 500 ml of test fluid required to immerse the spindles in a 600 ml beaker is not available. In practice, a smaller vessel may be used and the guard leg is removed. The Brookfield Viscometer/Rheometer will produce an

accurate and repeatable torque reading under any measurement circumstance. However, the conversion of this torque reading to centipoise will only be correct if the factor used was developed for those specific conditions. Brookfield has outlined a method for recalibrating a Brookfield Viscometer/Rheometer to any measurement circumstance in [More Solutions to Sticky Problems](#), Section 3.3.10. It is important to note that for many viscometer users the true viscosity is not as important as a repeatable day to day value. This repeatable value can be obtained without any special effort for any measurement circumstance. But, it should be known that this type of torque reading will not convert into a correct centipoise value when using a Brookfield factor if the boundary conditions are not those specified by Brookfield.

The guard leg is a part of the calibration check of the Brookfield LV and RV series Viscometer/Rheometer. Our customers should be aware of its existence, its purpose and the effect that it may have on data. With this knowledge, the viscometer user may make modifications to the recommended method of operation to suit their needs.

## APPENDIX F - VS-27Y Clamp Assembly

The VS-35Y Clamp assembly holds the DV-III+ on the upright rod and thus supports it on the base. The parts are shown in **Figure F1**.

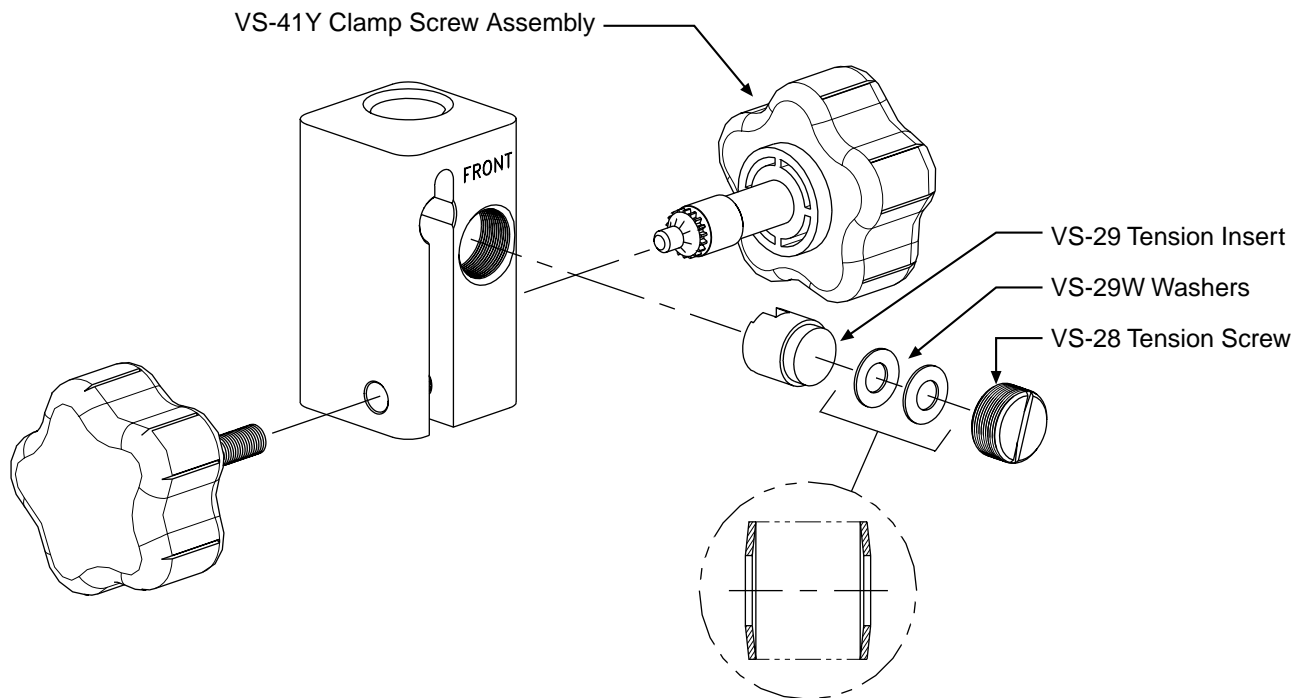

*Figure F1*

If the clamp is taken off the upright rod, the tension insert (Part No. VS-29) must be properly aligned for the clamp to fit back onto the upright rod.

When the tension insert (Part No. VS-29) is inserted, its slot must be in the vertical position parallel to the upright rod. If the slot is not in the correct position, the clamp will not slide down over the upright rod. Use a small screw driver or pencil to move it into the correct position. The VS-29W Belleville spring washers must face each other as illustrated. Adjust the VS-28 tension screw so that the clamp assembly is not loose on the upright rod.

## APPENDIX G - DV-III+ to Computer Command Set

The command set used to communicate with the DV-III+ is as follows:

Command	Format from Host	Response from Rheometer	Description
<b>E</b> (nable)	< <b>E</b> ><CR>	< <b>E</b> ><ss><CR>	Enable control circuitry
<b>R</b> (etrieve)	< <b>R</b> ><CR>	< <b>R</b> ><qqqq><tttt><ss><CR>	Retrieve data
<b>V</b> (elocity)	< <b>V</b> ><vvvvv><CR>	< <b>V</b> ><ss><CR>	Send speed
<b>I</b> (dentify)	< <b>I</b> ><CR>	< <b>I</b> ><dddd><mm><xx><CR>	Identify instrument
<b>Z</b> (ero)	< <b>Z</b> ><CR>	< <b>Z</b> ><zzzz><ss><CR>	Zero instrument
<b>Illegal String</b>	<????><CR>	<ss><CR>	Invalid command

Where:

- qqqq** = Transducer reading as **4** hex digits. The instrument should yield a reading of approximately **0400H** after zeroing at rest. Calculate % torque as follows:  

$$\% \text{ Torque} = \text{tval}/100$$
 where **tval** is the decimal equivalent of the received torque packet.
- zzzz** = Transducer reading representing the % torque zero offset as **4** hex digits. Calculate the decimal equivalent of the zero offset as follows:  

$$\text{Zero offset} = \text{tval}/100$$
 where **tval** is the decimal equivalent of the received zero offset packet. This value should now be subtracted from all future % torque readings retrieved using the **R** command.
- tttt** = Temperature reading as **4** hex digits. Calculate temperature in °C as follows:  

$$\text{Temperature } (^{\circ}\text{C}) = (\text{tval} - 4000)/40$$
 where **tval** is the decimal equivalent of the received temperature packet.
- vvvvv** = Stepper motor speed as **5** hex digits. The motor speed command from the Host is in units of **RPM**. To create a speed packet, multiply the desired decimal speed by 100 then convert the result to hexadecimal. All packets must be 5 characters long so for packets with less than that, pad the left side with zeros.  
**Example: To run the rheometer at 10 rpm**  

$$10 \text{ rpm} \times 100 = 1000 \text{ decimal} = 3\text{E8} \text{ hexadecimal};$$
**padding with zeros results in 003E8 so the command sent would be V003E8**
- ss** = Status Byte as **2** hex digits. (See Table 2)
- dddd** = The ASCII characters "**DV3+**".
- xx** = The firmware version of the instrument (i.e. version 4.1 is returned as **41**)
- mm** = These **2** characters will be used to represent the model for which the **DV-III+** is configured. The model decoding is shown in **Table 1** below.

Model	mm
LV	LV
2.5LV	4L
5LV	5L
RV	RV
RV	1R
RV	2R
HA	HA
2HA	3A
2.5HA	4A
HB	HB
2HB	3B
2.5HB	4B
5HB	5B

**Table 1 - DV-III+ Model Decoding**

Invalid cmd received	1	x	x	x	x	x	x	x
Diagnostics error	x	x	1	x	x	x	x	x
Motor circuit error	x	x	x	1	x	x	x	x
Auto-zero complete	x	x	x	x	1	x	x	x
Motor speed = 0	x	x	x	x	x	1	x	x
Motor On (Energized)	x	x	x	x	x	x	1	x
Control circuitry enabled	x	x	x	x	x	x	x	1

**Table 2 - Status Byte <ss> Definition**

**Note:** These values are not updated in the status byte when the listed condition occurs. They are made available when the computer next sends a command that includes the status byte in the response. The flags are cleared when the condition causing a flag to be set has been resolved or by re-enabling (E command) the **DV-III+**.

The **DV-III+** must first be issued the **E(nable)** command to enable control circuitry and ascertain its current status. The **DV-III+** will respond with an echo of the **E(nable)** command and will append the current status <ss> of the **DV-III+**. This 2-digit status byte will provide information as to the rheometer's internal working condition and capability to continue with or to accept new tasks.

The **Z(ero)** command is used to "zero" the **DV-III+** rheometer. The value returned <zzzz> is usually in the range of **03F0** hex to **0400** hex. This number should be retained and subtracted from every future returned torque reading to obtain the actual rheometer torque in percent.

The rheometer torque value and current temperature may be obtained by issuing the **R(etrieve)** command. The **DV-III+**'s response to this command is to echo the command: **R(etrieve)** and then follow it with **8** hex digits comprising the current values for the rheometer torque <qqqq> and the temperature probe reading <tttt>. Any control program written to exploit the output of the **DV-III+** will have to perform all the calculations required for viscosity, shear stress and shear rate, etc. The **DV-III+** will be directed to run at a given speed through the use of the **V(elocity)** command. Unlike the two previous commands, this command requires the parameter <vvvvv> which contains the value for the desired speed. The **DV-III+** responds to this command by repeating the command: **V(elocity)** and appending the status byte <ss> which will, amongst other things, inform the calling program as to whether the motor was turned on and/or whether the desired speed was attained.

## APPENDIX H - Fault Diagnosis and Troubleshooting

Listed are some of the more common problems that you may encounter while using your rheometer.

### ❑ Spindle Does Not Rotate

- ✓ Make sure the rheometer is plugged in.
- ✓ Check the voltage rating on your rheometer (115V or 220V); it must match the wall voltage.
- ✓ Make sure the motor is ON and the desired rpm is selected.

### ❑ Spindle Wobbles When Rotating or Looks Bent

- ✓ Make sure the spindle is tightened securely to the rheometer coupling.
- ✓ Check the straightness of all other spindles; replace if bent.
- ✓ Inspect rheometer coupling and spindle coupling mating areas and threads for dirt; clean threads on spindle coupling with a 3/56 left-hand tap.
- ✓ Inspect rheometer coupling threads for wear; if the threads are worn, the unit needs service (see Appendix I). Check to see if spindles rotate eccentrically or wobble. There is an allowable runout of 1/32-inch in each direction (1/16-inch total) when measured from the bottom of the spindle rotating in air.
- ✓ Check to see if the rheometer coupling appears bent; if so, the unit is in need of service (see Appendix I, “How to Return Your Rheometer”).

If you continue to experience problems with your rheometer, follow this troubleshooting section to help isolate potential problems.

### ❑ Perform an Oscillation Check

This check verifies the mechanical condition of the shaft and bearing assembly in the rheometer. Auto zero the instrument according to directions in the manual. With the spindle removed, motor off and the display showing % torque perform the oscillation check as follows:

- 1) Manually lift the spindle nut and rotate counter clockwise until the display shows 20%-30% deflection.
- 2) Release the spindle nut.
- 3) The % display should return to 0.0%  $\pm$ 0.1%.

### ❑ **Viscometer Will Not Return to Zero**

- ✓ Rheometer is not level
  - Adjust the laboratory stand
- ✓ Pivot point or jewel bearing faulty
  - Perform calibration check
  - Contact Brookfield Engineering Laboratories, Inc. or your Brookfield dealer for repair.

### ❑ **Inaccurate Readings**

- ✓ Verify spindle, speed and model selection.
- ✓ If % readings are under-range (less than 10%), the units display (% , cP, D/cm<sup>2</sup>, 1/sec) will flash; change spindle and/or speed.
- ✓ “EEEE” on the digital display means the unit is over-range (greater than 100%); reduce speed and/or change spindle.
- ✓ Verify test parameters: temperature, container, volume, method. Refer to:
  - “More Solutions to Sticky Problems”, Section III
  - DV-III+ Digital Rheometer Operating Instructions, Appendix C, “Variables in Viscosity Measurements.”
- ✓ Perform a calibration check; follow the instructions in Appendix E.
- ✓ Verify tolerances are calculated correctly.
- ✓ Verify the calibration check procedures were followed exactly.

If the unit is found to be out of tolerance, the unit may be in need of service. See Appendix I for details on “How to Return Your Rheometer.”

### ❑ **Display Reading Will Not Stabilize**

- ✓ Special characteristic of sample fluid. There is no problem with the rheometer.
  - Refer to Appendix C
- ✓ Check for erratic spindle rotation
  - Verify power supply
  - Contact Brookfield Engineering Laboratories, Inc. or your Brookfield dealer for repair.
- ✓ Bent spindle or spindle coupling.
  - Contact Brookfield Engineering Laboratories, Inc. or your Brookfield dealer for repair.
- ✓ Temperature fluctuation in sample fluid
  - Use temperature bath control.

### **No Recorder Response**

- ✓ Be sure the rheometer is not at ZERO reading.
- ✓ Be sure the recorder is ON and not on STANDBY.
- ✓ Verify the range settings.
- ✓ Check cable leads for clean connection.

### **Recorder Pen Moves in Wrong Direction**

- ✓ Output polarity reversed
  - Reverse leads

### **Rheometer Will Not Communicate with PC**

- ✓ Contact Brookfield or an authorized dealer with the following information:
  - When the communication error occurs
  - The exact text displayed when the error occurs
  - Computer hardware details including processor speed, RAM, network cards, modems, etc.
  - Computer operating system
  - Total number of COM ports
  - List of Brookfield equipment attached to COM ports and their corresponding COM port number

## APPENDIX I - Warranty Repair and Service

### Warranty

Brookfield Viscometers are guaranteed for one year from date of purchase against defects in materials and workmanship. They are certified against primary viscosity standards traceable to the National Institute of Standards and Technology (NIST). The Viscometer must be returned to **Brookfield Engineering Laboratories, Inc.** or the Brookfield dealer from whom it was purchased for no charge warranty service. Transportation is at the purchaser's expense. The Viscometer should be shipped in its carrying case together with all spindles originally provided with the instrument.

For repair or service **in the United States**, return to:

**Brookfield Engineering Laboratories, Inc.**

11 Commerce Boulevard  
Middleboro, MA 02346 U.S.A.

Telephone: (508) 946-6200 FAX: (508) 946-6262  
e-mail: sales@brookfieldengineering.com website: www.brookfieldengineering.com

For repair or service **outside the United States**, consult **Brookfield Engineering Laboratories, Inc.** or the dealer from whom you purchased the instrument.

For repair or service in the **United Kingdom**, return to:

**Brookfield Viscometers Limited**

1 Whitehall Estate  
Flex Meadow  
Pinnacles West  
Harlow, Essex CM19 5TJ, United Kingdom

Telephone: (44) 27/945 1774 FAX: (44) 27/945 1775  
e-mail: sales@brookfield.co.uk website: www.brookfield.co.uk



For repair or service in **Germany**, return to:

**Brookfield Engineering Laboratories Vertriebs GmbH**

Barbarossastrasse 3  
D-73547 Lorch, Germany

Telephone: (49) 7172/927100 FAX: (49) 7172/927105  
e-mail: info@brookfield-gmbh.de

## Packaging Instructions to Return a Viscometer for Repair or Calibration

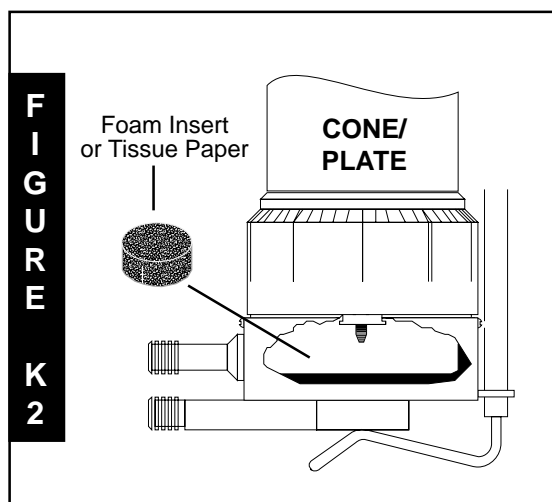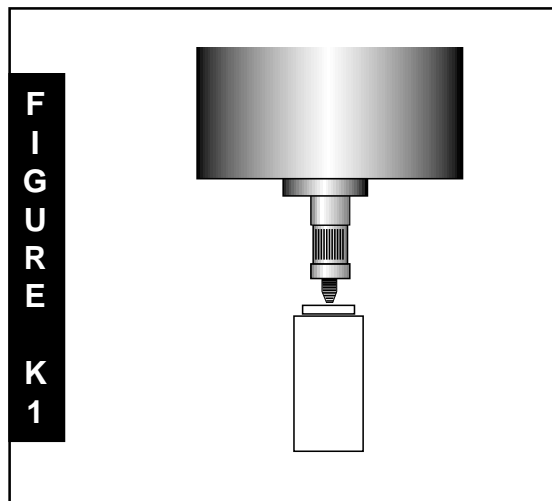

- Remove and return all spindles (properly packed for shipping).
- Clean excess testing material off the instrument.
- Include MSDS sheets for all materials tested with this instrument.
- Support pointer shaft with white, nylon shipping cap, as shown in Figure K1, or with white plastic shipping cap originally supplied with instrument.
- Pack the instrument in its original case. Cases are available for immediate shipment from Brookfield. If the case is not available, take care to wrap the instrument with enough material to support it. Avoid using foam peanuts or shredded paper.
- DO NOT send the laboratory stand unless there is a problem with the upright rod, clamp or base. If there is a problem with the stand, remove the upright rod from the base and individually wrap each item to avoid contact with the instrument. Do not put lab stand in viscometer carrying case.
- Fill out a copy of the Viscometer Information Sheet (on the following page) with as much information as possible to help expedite your service *or* include a memo indicating the type of problem you are experiencing or the service you need performed. Please also include a purchase order number for us to bill against.
- Package the instrument and related items in a strong box for shipping. Mark the outside of the box with handling instructions.

**Example:** "Handle with Care" or  
"Fragile - Delicate Instrument"

For cone/plate instruments, please remove the cone spindle and carefully pack in place in the shipping case. If available, use the original foam insert or roll up one sheet of tissue paper (or similar) and place between the spindle coupling and cup assembly (see Figure K2). This will help prevent damage in shipping.

Providing us with the following information will help us to service your equipment more quickly and efficiently. Please photocopy, fill out and return a copy of this form with your instrument.

Brookfield recommends that all viscometers be returned for annual calibration to ensure that your equipment continues to provide the same accuracy you have come to expect from Brookfield products.

<b>1</b>	<b>VISCOMETER INFORMATION</b>
	<b>Serial Number:</b> _____ <b>Date:</b> _____ <b>Model:</b> _____

<b>2</b>	<b>COMPANY INFORMATION</b>
	Company: _____ Primary User: _____ Telephone: _____ Fax: _____ P.O. Number: _____ (to cover repair and shipping) Billing Address: _____ Shipping Address: _____ _____ _____ Return Shipment Instructions: UPS ground <input type="checkbox"/> UPS Next Day <input type="checkbox"/> UPS 2nd Day <input type="checkbox"/> Federal Express <input type="checkbox"/> _____ (Federal Express Account Number required)

<b>3</b>	<b>SERVICE INFORMATION</b>
	Operating Conditions (Spindle; Speed; Viscosity Range; Temp. Control; Temperature of Sample): _____ _____ Did you contacting Brookfield before returning this instrument? Y <input type="checkbox"/> N <input type="checkbox"/> If yes, whom did you contact? _____ Description/Symptoms of Present Problem/Malfunction (please list all): _____ _____ Time Since Last Serviced (if known): _____ Before & After Calibration Check? Y <input type="checkbox"/> N <input type="checkbox"/> Other Comments: _____ _____ _____

<b>STEPS:</b> <input type="checkbox"/> Return the Viscometer to the attention of the Repair Department at the address above. <input type="checkbox"/> Package the Viscometer for shipment, as outlined on the previous page. <input type="checkbox"/> Include a purchase order or purchase order number with this form.
---

*This tear-off sheet is a typical example of recorded test data. Please photocopy and retain this template so that additional copies may be made as needed.*

<h1 style="margin: 0;">VISCOSITY TEST REPORT</h1>										FOR:
TEST INFORMATION:										
DATE: _____ BY: _____										
MODEL: _____										
SPINDLE: _____										
RPM: _____										
DIAL READING % TORQUE: _____										
FACTOR: _____										
VISCOSITY cP: _____										
SHEAR RATE: _____										
TEMP °C: _____										
TIME: _____										
NOTES:										
CONCLUSIONS:										

BROOKFIELD ENGINEERING LABORATORIES, INC. • 11 Commerce Boulevard • Middleboro, MA 02346 • TEL: 508-946-6200 or 800-628-8139 (USA ex MA) FAX: 508-946-6262 • www.brookfieldengineering.com